

Protein significance analysis of mass spectrometry-based proteomics experiments with R and MSstats (v3.18.1) or later

Meena Choi & Tsung-Heng Tsai

February 26, 2020

Contents

1. Statistical relative protein quantification: SRM, DDA and DIA experiments	2
1.1 Applicability	2
1.2 Statistical functionalities	2
1.3 Interoperability with existing computational tools	2
1.4 Availability	3
1.5 Overview of the functionalities	4
1.6 Troubleshooting	5
2. Allowable data formats	5
2.1 SRM with stable isotope labeled reference peptides	5
2.2 Label-free DDA	9
2.2 Label-free DIA	10
3. Prerequisites and setting for MSstats analysis	10
4. DDA analysis with MSstats	12
4.1 General workflow for DDA	12
4.2 Suggested workflow with Skyline output for DDA	33
4.3 Suggested workflow with MaxQuant output for DDA	37
4.4 Suggested workflow with Progenesis output for DDA	39
4.5 Suggested workflow with Proteome Discoverer output for DDA	44
4.6 Suggested workflow with OpenMS output for DDA	49
5. DIA analysis with MSstats	51
5.1 Suggested workflow with Skyline output for DIA	51
5.2 Suggested workflow with Spectronaut output for DIA	51
5.3 Suggested workflow with DIA-Umpire output for DIA	55
5.4 Suggested workflow with OpenSWATH output for SWATH	59
6. SRM analysis with MSstats	65
6.1 Suggested workflow for SRM	65
Reference	68

1. Statistical relative protein quantification: SRM, DDA and DIA experiments

MSstats is an open-source R-based package for statistical relative quantification of peptides and proteins in mass spectrometry-based proteomic experiments. This document describes **MSstats**, the most recent version of the package, and its use through the command line.

1.1 Applicability

MSstats version 3.0 and above is applicable to multiple types of sample preparation, including label-free workflows, workflows that use stable isotope labeled reference proteins and peptides, and workflows that use fractionation. It is applicable to targeted Selected Reaction Monitoring (SRM), Data-Dependent Acquisition (DDA or shotgun), and Data-Independent Acquisition (DIA or SWATH-MS). It is applicable to experiments that make arbitrary complex comparisons of experimental conditions or times.

MSstats is not applicable to experiments that compare multiple metabolically labeled endogenous samples within a same run, such as experiments with iTRAQ labeling or TMT labeling. These experiments are supported by **MSstatsTMT**, which is a sibling package. Please check **MSstatsTMT** in Bioconductor.

1.2 Statistical functionalities

MSstats version 3.0 and above performs three analysis steps. The first step, *data processing, visualization, and run-level summarization*, transforms, normalizes and summarizes the intensities of the peaks per MS run and per protein, and generates workflow-specific and customizable numeric summaries for data visualization and quality control.

The second step, *statistical modeling and inference*, automatically detects the experimental design (e.g. group comparison, paired design or time course, presence of labeled reference peptides or proteins) from the data. It then reflects the experimental design and the type of spectral acquisition strategy, and fits an appropriate linear mixed model by means of **lm** and **lmer** functionalities in R. The model is used to detect differentially abundant proteins or peptides, or to summarize the protein or peptide abundance in a single biological replicate or condition (that can be used, e.g. as input to clustering or classification).

The third step, *statistical experimental design*, views the dataset being analyzed as a pilot study of a future experiment, utilizes the variance components of the current dataset, and calculates the minimal number of replicates necessary in the future experiment to achieve a pre-specified statistical power.

1.3 Interoperability with existing computational tools

MSstats takes as input data in a tabular .csv format, which can be generated by any spectral processing tool such as Skyline (MacLean et al. 2010), MaxQuant (Cox and Mann 2008), Progenesis QI (Nonlinear dynamics/Waters), Proteome Discoverer (Thermo Scientific), MultiQuant (Applied Biosystems), OpenMS (Sturm et al. 2008), SuperHirn (Mueller et al. 2007), OpenSWATH (Röst et al. 2014), Spectronaut (Biognosys), or DIA-Umpire (Tsou et al. 2015). The functions to convert the required format from several processing tools are available from **MSstats** v3.6. Details are in the section below.

For statistics experts, **MSstats** 3.0 and above satisfies the interoperability requirements of Bioconductor. The command line-based workflow is partitioned into a series of independent steps, that facilitate the development and testing of alternative statistical approaches. It complies with the maintenance and documentation requirements of Bioconductor.


MSstats 3.0 and above is available as an external tool within Skyline. The external tool support within Skyline manages **MSstats** installation, point-and-click execution, parameter collection in Windows forms and

output display. Skyline manages the annotations of the experimental design, and the processing of raw data. It outputs a custom report, that is fed as a single stream input into **MSstats**. This design buffers proteomics users from the details of the R implementation, while enabling rigorous statistical modeling. Also, **MSstat** can be combined with an OpenMS preprocessing pipeline (e.g. in KNIME). The OpenMS experimental design is used to present the data in an MSstats-conformant way for the analysis. Details are available in OpenMS tutorial.

1.4 Availability

MSstats is available under the Artistic-2.0 license at msstats.org. **MSstats** as an external tool for Skyline is available at <http://proteome.gs.washington.edu/software/Skyline/tools.html>. **MSstats** is now also available in Bioconductor. The most recent version of the package is available at msstats.org or MSstats GitHub. We suggest to use that if possible. The versioning of the main package is updated several times a year, to synchronise with the Bioconductor release.

1.5 Overview of the functionalities

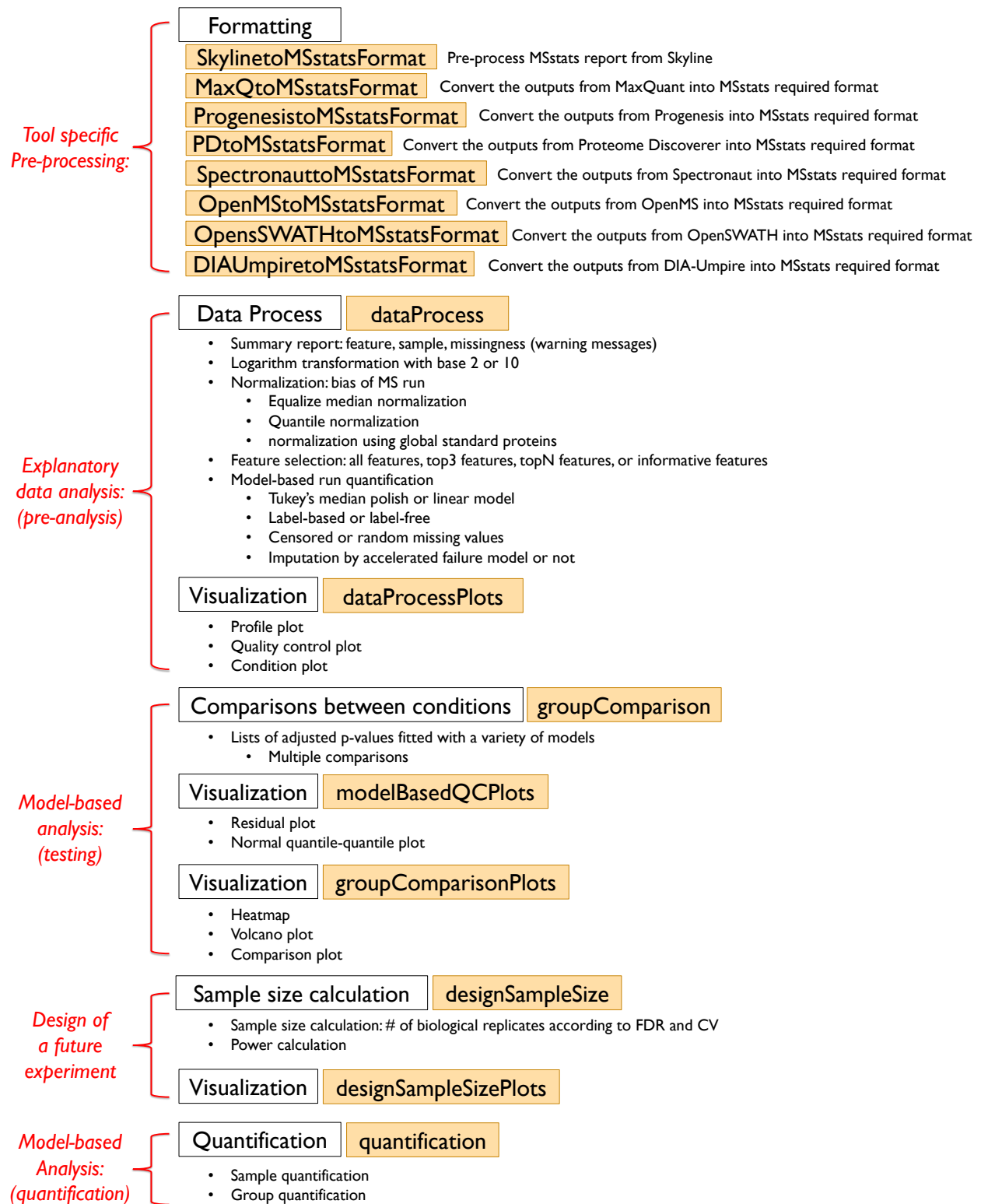

1.6 Troubleshooting

To help troubleshoot potential problems with installation or functionalities of **MSstats**, a progress report is generated in a separate log file, *msstats.log*. The file includes information on the R session (R version, loaded software libraries), options selected by the user, checks of successful completion of intermediate analysis steps, and warning messages. If the analysis produces an error, the file contains suggestions for possible reasons for the errors. If a file with this name already exists in working directory, a suffix with a number will be appended to the file name. In this way a record of all the analyses is kept. Please see the file **KnownIssues-Skyline-MSstatsV3.6.pdf** on the “Installation” section of “MSstats” page in msstats.org for a list of known issues and possible solutions for installation problem of MSstats external tool in Skyline

2. Allowable data formats

2.1 SRM with stable isotope labeled reference peptides

2.1.1 10-column format

MSstats performs statistical analysis steps, that follow peak identification and quantitation. Therefore, input to **MSstats** is the output of other software tools (such as Skyline or MultiQuant) that read raw spectral files and identify and quantify spectral peaks. The preferred structure of data for use in **MSstats** is a .csv file in a “long” format with 10 columns representing the following variables: **ProteinName**, **PeptideSequence**, **PrecursorCharge**, **FragmentIon**, **ProductCharge**, **IsotopeLabelType**, **Condition**, **BioReplicate**, **Run**, **Intensity**. The variable names are fixed, but are case-insensitive.

- (a) **ProteinName**: This column needs information about Protein id. Statistical analysis will be done separately for each unique label in this column. For peptide-level modeling and analysis, use peptide id in this column.

(b)-(e) **PeptideSequence**, **PrecursorCharge**, **FragmentIon**, **ProductCharge**: The combination of these 4 columns defines a *feature* of a protein (in SRM experiments, it is a transition that is identified and quantified across runs). If the information for one or several of these columns is not available, please do not discard these columns but use a single fixed value across the entire dataset. For example, if the original raw data does not contain the information of **ProductCharge**, assign the value 0 to the entries in the column **ProductCharge** for the entire dataset. If the peptide sequences should be distinguished based on post-translational modifications, this column can be renamed to **PeptideModifiedSequence**. For example, this allows us to use the **PeptideModifiedSequence** column from the Skyline report.

- (f) **IsotopeLabelType**: This column indicates whether this measurement is based on the endogenous peptides (use “L”) or labeled reference peptides (use “H”).
- (g) **Condition**: For group comparison experiments, this column indicates groups of interest (such as “Disease” or “Control”). For time-course experiments, this column indicates time points (such as “T1”, “T2”, etc). If the experimental design contains both distinct groups of subjects and multiple time points per subject, this column should indicate a combination of these values (such as “Disease_T1”, “Disease_T2”, “Control_T1”, “Control_T2”, etc.).
- (h) **BioReplicate**: This column should contain a unique identifier for each biological replicate in the experiment. For example, in a clinical proteomic investigation this should be a unique patient id. Patients from distinct groups should have distinct ids. **MSstats** does not require the presence of technical replicates in the experiment. If the technical replicates are present, all samples or runs from a same biological replicate should have a same id. **MSstats** automatically detects the presence of technical replicates and accounts for them in the model-based analysis.

- (i) **Run:** This column contains the identifier of a mass spectrometry run. Each mass spectrometry run should have a unique identifier, regardless of the origin of the biological sample. In SRM experiments, if all the transitions of a biological or a technical replicate are split into multiple “methods” due to the technical limitations, each method should have a separate identifier. When processed by Skyline, distinct values of runs correspond to distinct input file names. It is possible to use the actual input file names as values in the column **Run**.
- (j) **Intensity:** This column should contain the quantified signal of a feature in a run without any transformation (in particular, no logarithm transform). The signals can be quantified as the peak height or the peak of area under curve. Any other quantitative representation of abundance can also be used.

An example of an acceptable input dataset is shown below. This example dataset is from an SRM experiment with stable isotope labeled reference peptides. The dataset is stored in a .csv file in a “long” format. Each row corresponds to a single intensity. More details on assigning the values of **Condition**, **BioReplicate** and **Run**, depending on the structure of the experimental design, are given below.

	A	B	C	D	E	F	G	H	I	J
1	ProteinName	PeptideSequence	PrecursorCharge	FragmentIon	ProductCharge	IsotopeLabelType	Condition	BioReplicate	Run	Intensity
2	ACEA	EILGHEIFFDWELP	3	y3	0	H	1	ReplA	1	66472.3847
3	ACEA	EILGHEIFFDWELP	3	y3	0	L	1	ReplA	1	5764.16228
4	ACEA	EILGHEIFFDWELP	3	y4	0	H	1	ReplA	1	101005.166
5	ACEA	EILGHEIFFDWELP	3	y4	0	L	1	ReplA	1	61.65238
6	ACEA	EILGHEIFFDWELP	3	y5	0	H	1	ReplA	1	90055.4993
7	ACEA	EILGHEIFFDWELP	3	y5	0	L	1	ReplA	1	472.691803
8	ACEA	TDSEAATLISSTID	2	y10	0	H	1	ReplA	1	43506.5425
9	ACEA	TDSEAATLISSTID	2	y10	0	L	1	ReplA	1	217.203553
10	ACEA	TDSEAATLISSTID	2	y7	0	H	1	ReplA	1	68023.0377
11	ACEA	TDSEAATLISSTID	2	y7	0	L	1	ReplA	1	725.284308
12	ACEA	TDSEAATLISSTID	2	y8	0	H	1	ReplA	1	68276.0489
13	ACEA	TDSEAATLISSTID	2	y8	0	L	1	ReplA	1	243.658527

2.1.2 Assigning the values of Condition, BioReplicate and Run

The values of **Condition**, **BioReplicate**, **Run** depend on the design of the specific experiment.

1) Group comparison In a group comparison design, the conditions (e.g., disease states) are profiled across **non-overlapping sets of biological replicates** (i.e., subjects). In this example there are 2 conditions, Disease and Control (in general the number of conditions can vary). There are 3 subjects (i.e., biological replicates) per condition (in general an equal number of replicates per condition is not required). Each subject has 2 technical replicate runs (in general technical replicates are not required, and their number per sample may vary). Overall, in this example there are $2 \times 3 \times 2 = 12$ mass spectrometry runs.

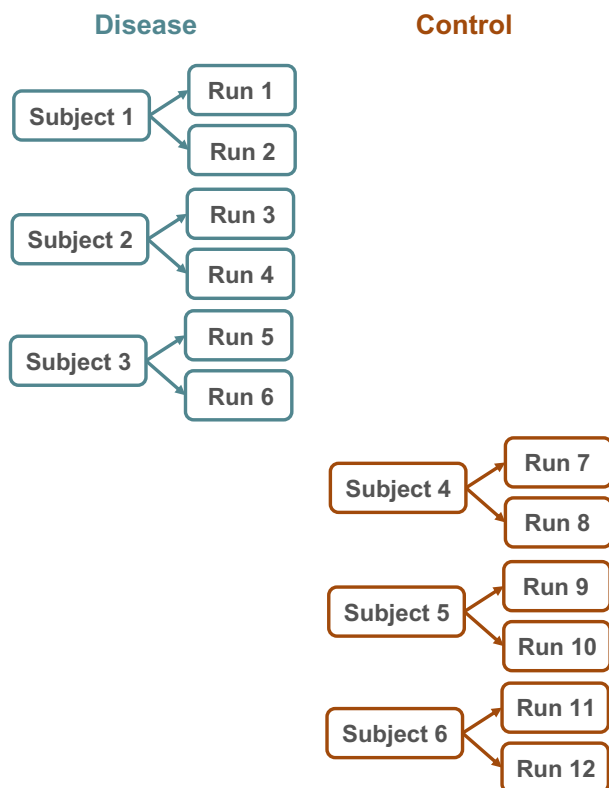


Table below shows the values of the columns **Condition**, **BioReplicate** and **Run** for this situation. It is important to note two things. First, the order of subjects and conditions in the experiment should be randomized, and run id does not need to represent the order of spectral acquisition. Second, the values of the columns are repeated for every quantified transition. For example, if in each run the experiment quantifies 50 endogenous transitions and 50 labeled reference counterparts, then the input file has $12 \times 50 \times 2 = 1200$ lines. When a feature intensity is missing in a run, the data structure should contain a separate row for each missing value. The rows should include all the information (from **ProteinName** to **Run**), and indicate missing intensities with NA.

Condition	BioReplicate	Run
Disease	Subject1	1
Disease	Subject1	2
Disease	Subject2	3
Disease	Subject2	4
Disease	Subject3	5
Disease	Subject3	6
Control	Subject4	7
Control	Subject4	8
Control	Subject5	9
Control	Subject5	10
Control	Subject6	11
Control	Subject6	12

2) Time course The important feature of a time course experimental design is that **a same subject (i.e., biological replicate) is repetitively measured across multiple time points**. In this example there are 2 time points, Time1 and Time2 (in general the number of times can vary). There are 4 subjects (i.e., biological replicates) measured across times (in general an equal number of times per replicate is not

required). There are no technical replicates (in general the number of technical replicates per sample may vary). Overall, in this example there are $2 \times 4 \times 1 = 8$ mass spectrometry runs.

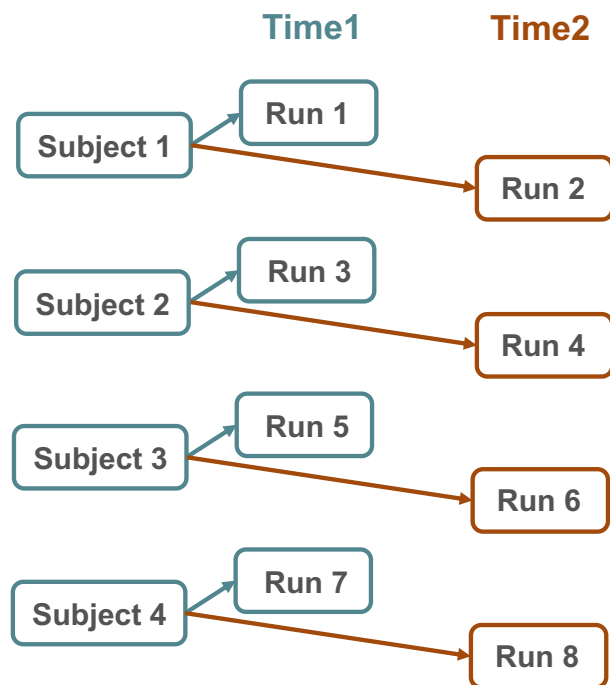


Table below shows the values of the columns **Condition**, **BioReplicate** and **Run** for this situation. Comments on the order of the runs, on the number of lines in the input data structure, and on the handling of missing peak intensities are as in the group comparison design.

Condition	BioReplicate	Run
Time1	Subject1	1
Time2	Subject1	2
Time1	Subject2	3
Time2	Subject2	4
Time1	Subject3	5
Time2	Subject3	6
Time1	Subject4	7
Time2	Subject4	8

3) Paired design Another frequently used experimental design is a *paired design*, where measurements from multiple conditions (such as healthy biopsy and disease biopsy) are taken from a same subject. The statistical model for this experimental design is the same as in the time course experiment, however the values in the columns of the input data may have a different appearance. In this example there are 2 subjects, PatientA and PatientB (in general the number of patients can vary). There are two conditions per subject, BiopsyHealthy and BiopsyTumor (in general the number of conditions per subject can exceed two). In this example there are 3 technical replicates of each type (in this example, the technical replicates are biopsies; in general these can also be replicate sample preparations or replicate mass spectrometry runs). Overall, in this example there are $2 \times 2 \times 3 = 12$ mass spectrometry runs.

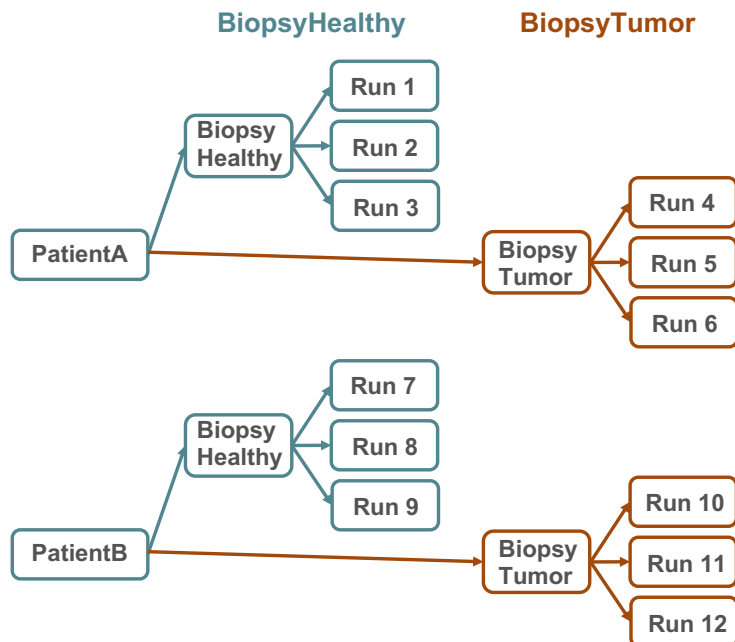


Table below shows the values of the columns **Condition**, **BioReplicate** and **Run** for this situation. Comments on the order of the runs, on the number of lines in the input data structure, and on the handling of missing peak intensities are as in the group comparison design.

Condition	BioReplicate	Run
BiopsyHealthy	PatientA	1
BiopsyHealthy	PatientA	2
BiopsyHealthy	PatientA	3
BiopsyTumor	PatientA	4
BiopsyTumor	PatientA	5
BiopsyTumor	PatientA	6
BiopsyHealthy	PatientB	7
BiopsyHealthy	PatientB	8
BiopsyHealthy	PatientB	9
BiopsyTumor	PatientB	10
BiopsyTumor	PatientB	11
BiopsyTumor	PatientB	12

2.2 Label-free DDA

For label-free DDA experiments the required input is the 10-column format, the same as described in section 2.1 for SRM experiments. In DDA experiments spectral features are defined as peptide ions, which are identified and quantified across runs. Since for label-free DDA experiments some of the columns **PeptideSequence**, **PrecursorCharge**, **FragmentIon**, and **ProductCharge** are not relevant, these columns will have a constant fixed value (such as **NA**) across the entire dataset. Furthermore, the column **IsotopeLabelType** will be set to “L” for the entire dataset.

ProteinName	PeptideSequence	PrecursorCharge	Fragmention	ProductCharge	IsotopeLabelType	Condition	BioReplicate	Run	Intensity
bovine	S.PVDIDTK	5	NA	NA	L	C1	1	1	2636791.5
bovine	S.PVDIDTK	5	NA	NA	L	C1	1	2	1992418.5
bovine	S.PVDIDTK	5	NA	NA	L	C1	1	3	1982146.38
bovine	S.PVDIDTK	5	NA	NA	L	C2	1	4	5019594
bovine	S.PVDIDTK	5	NA	NA	L	C2	1	5	4560467.5
bovine	S.PVDIDTK	5	NA	NA	L	C2	1	6	3627848.75
bovine	S.PVDIDTK	5	NA	NA	L	C5	1	13	145511.83
bovine	S.PVDIDTK	5	NA	NA	L	C5	1	14	291829.69
bovine	S.PVDIDTK	5	NA	NA	L	C6	1	16	786667.38
bovine	S.PVDIDTK	5	NA	NA	L	C6	1	17	705295.31
bovine	S.PVDIDTK	5	NA	NA	L	C6	1	18	453448.78
bovine	S.PVDIDTK	5	NA	NA	L	C3	1	7	NA

2.2 Label-free DIA

For label-free DIA experiments, the required input is the 10-column format, the same as described in section 2.1 for SRM experiments. The values of the required columns can be extracted from the output of signal processing software such as Skyline or OpenSWATH. By default, the combination of the values in the columns `PeptideSequence`, `PrecursorCharge`, `Fragmention`, `ProductCharge` uniquely identifies each spectral feature (i.e., a fragment ion identified and quantified across multiple runs). If the signal processing software does not provide the information on some of these columns but provides a unique feature identifier, it is possible to use this unique identifier instead of one of these columns. Furthermore, the column `IsotopeLabelType` is set to “L” for the entire dataset.

An example dataset is shown below. In this example, feature id generated by OpenSWATH is used instead of `ProductCharge` to uniquely characterize each feature.

ProteinName	PeptideSequence	PrecursorCharge	Fragmention	ProductCharge	IsotopeLabelType	Condition	BioReplicate	Run	Intensity
350748	TPPAAVLLK	2	y7	109401	L	2	1	3	257486
350748	TPPAAVLLK	2	y7	109401	L	2	2	4	141159
350748	TPPAAVLLK	2	y7	109401	L	1	1	1	452908
350748	TPPAAVLLK	2	y7	109401	L	1	2	2	348222
515084	NIC[160]VNAIAPGFIESDMTGVLPEK	3	y3	7717	L	2	1	3	12753
515084	NIC[160]VNAIAPGFIESDMTGVLPEK	3	y3	7717	L	2	2	4	12857
515084	NIC[160]VNAIAPGFIESDMTGVLPEK	3	y3	7717	L	1	1	1	89652
515084	NIC[160]VNAIAPGFIESDMTGVLPEK	3	y3	7717	L	1	2	2	76724
515084	MVNEAIESLGSIDVLVNNAGITNDK	3	y9	57971	L	2	1	3	2052
515084	MVNEAIESLGSIDVLVNNAGITNDK	3	y9	57971	L	2	2	4	1050
515084	MVNEAIESLGSIDVLVNNAGITNDK	3	y9	57971	L	1	1	1	10772
515084	MVNEAIESLGSIDVLVNNAGITNDK	3	y9	57971	L	1	2	2	10516

3. Prerequisites and setting for MSstats analysis

MSstats is an R-based package. It is assumed that you already have R installed. You can install MSstats from Bioconductor:

```
if (!requireNamespace("BiocManager", quietly = TRUE))
  install.packages("BiocManager")
```

```
BiocManager::install("MSstats")
```

Once you have the package installed, load MSstats into an R session and verify that you have the correct version. Note that in order to use MSstats, the package needs to be loaded every time you restart R.

```
library('MSstats', warn.conflicts = F, quietly = T, verbose = F)
?MSstats
sessionInfo()
```

```
## R version 3.6.2 (2019-12-12)
## Platform: x86_64-apple-darwin15.6.0 (64-bit)
## Running under: macOS Catalina 10.15.3
##
## Matrix products: default
## BLAS:   /Library/Frameworks/R.framework/Versions/3.6/Resources/lib/libRblas.0.dylib
## LAPACK: /Library/Frameworks/R.framework/Versions/3.6/Resources/lib/libRlapack.dylib
##
## locale:
## [1] en_US.UTF-8/en_US.UTF-8/en_US.UTF-8/C/en_US.UTF-8/en_US.UTF-8
##
## attached base packages:
## [1] stats      graphics  grDevices  utils      datasets  methods    base
##
## other attached packages:
## [1] MSstats_3.19.5
##
## loaded via a namespace (and not attached):
## [1] gtools_3.8.1      statmod_1.4.34    minpack.lm_1.2-1
## [4] tidyselect_1.0.0  xfun_0.12         reshape2_1.4.3
## [7] purrr_0.3.3       splines_3.6.2     lattice_0.20-40
## [10] colorspace_1.4-1  vctrs_0.2.3       generics_0.0.2
## [13] doSNOW_1.0.18     htmltools_0.4.0   snow_0.4-3
## [16] yaml_2.2.1        marray_1.64.0     survival_3.1-8
## [19] rlang_0.4.4       pillar_1.4.3      nloptr_1.2.1
## [22] glue_1.3.1        plyr_1.8.5        foreach_1.4.8
## [25] lifecycle_0.1.0   stringr_1.4.0     munsell_0.5.0
## [28] gtable_0.3.0      caTools_1.18.0    codetools_0.2-16
## [31] evaluate_0.14     knitr_1.28        parallel_3.6.2
## [34] preprocessCore_1.48.0 broom_0.5.4       Rcpp_1.0.3
## [37] KernSmooth_2.23-16 scales_1.1.0      backports_1.1.5
## [40] gdata_2.18.0      limma_3.42.0     lme4_1.1-21
## [43] ggplots_3.0.3     ggplot2_3.2.1     digest_0.6.25
## [46] stringi_1.4.6     ggrepel_0.8.1     dplyr_0.8.4
## [49] grid_3.6.2        bitops_1.0-6      tools_3.6.2
## [52] magrittr_1.5       lazyeval_0.2.2    tibble_2.1.3
## [55] crayon_1.3.4      tidyr_1.0.2       pkgconfig_2.0.3
## [58] MASS_7.3-51.5     Matrix_1.2-18     data.table_1.12.8
## [61] assertthat_0.2.1  minqa_1.2.4       rmarkdown_2.1
## [64] iterators_1.0.12  R6_2.4.1          boot_1.3-24
## [67] nlme_3.1-144      compiler_3.6.2
```

Finally, set the working directory to where you saved files. Note that you may have a different path on your computer from the example.

```
setwd('/Users/meenachoi/Dropbox/MSstats_GitHub_document/MSstats_v3.18.1')
```

You can check your working directory by:

```
getwd()
```

```
## [1] "/Users/meenachoi/Dropbox/MSstats_GitHub_document/MSstats_v3.18.1"
```

4. DDA analysis with MSstats

4.1 General workflow for DDA

This section describes a typical workflow for DDA analysis with **MSstats**. Controlled mixture DDA data will be used for demonstration. This dataset is available as an example data (**DDARawData**) in **MSstats**. Also the csv file for the same dataset, **RawData.DDA.csv**, is available in **MSstats** material GitHub in the folder named ‘example dataset/DDA_controlledMixture2009’. It is processed by Superhirn. (original reference link)

4.1.1 Preparing the data for MSstats input

The first step in using the **MSstats** is to format the data as described in Section 2. **DDARawData** is already formatted for **MSstats** input.

```
# Check the first 6 rows in DDARawData
head(DDARawData)
```

```
## ProteinName PeptideSequence PrecursorCharge FragmentIon ProductCharge
## 1 bovine S.PVDIDTK_5 5 NA NA
## 2 bovine S.PVDIDTK_5 5 NA NA
## 3 bovine S.PVDIDTK_5 5 NA NA
## 4 bovine S.PVDIDTK_5 5 NA NA
## 5 bovine S.PVDIDTK_5 5 NA NA
## 6 bovine S.PVDIDTK_5 5 NA NA
## IsotopeLabelType Condition BioReplicate Run Intensity
## 1 L C1 1 1 2636792
## 2 L C1 1 2 1992418
## 3 L C1 1 3 1982146
## 4 L C2 1 4 5019594
## 5 L C2 1 5 4560468
## 6 L C2 1 6 3627849
```

4.1.2 Processing the data

Normalizing and summarizing data with dataProcess After reading the datasets, **MSstats** performs 1) logarithm transformation of **Intensity** column, 2) normalization, 3) feature selection, (all features vs subset of features), 4) imputation for censored missing value, which are below the cutoff and undetectable, 5) run-level summarization.

To get started with this function, visit the help section of **dataProcess** first:

```
?dataProcess
```

NOTE At the logarithm transformation step, zero value in **Intensity** is problematic. When **Intensity**=0, **Inf** is the output from logarithm transformed intensities. Also, logarithm transformed intensities, when **Intensity** < 1, are negative values and it can make overestimated between log fold change. Therefore, logarithm transformed intensities for original intensity between 0 and 1 will be replaced with zero value after normalization.

Default normalization and summarization options **dataProcess** provides a variety of options in consideration of different experimental protocols. Default values for all options are our suggestion for general cases. However, the default options may not be appropriate for all possible scenarios. It is important to understand their underlying assumption to avoid misuse. Below is the additional explanation for main options.

- **logTrans** : logarithm transformation with base 2 (default) of **Intensity** column.
- **Normalization** :
 - ‘**equalizeMedians**’ : The default option for normalization is **equalizeMedians**, where all the intensities in a run are shifted by a constant, to equalize the median of intensities across runs for label-free experiment. This normalization method is appropriate when we can assume that the majority of proteins do not change across runs. *Be cautious when using the **equalizeMedians** option for a label-free DDA dataset with only a small number of proteins.* For label based experiment, **equalizeMedians** equalizes the median of reference intensities across runs and is generally proper even for a dataset with a small number of proteins.
 - ‘**globalStandards**’ : Instead, if you have a spiked in standard, you may set this to **globalStandards** and define the standard with **nameStandards** option.
 - ‘**quantile**’ : The distribution of all the intensities in each run will become the same across runs for label-free experiment. For label-based experiment, the distribution of all the reference intensities will be become the same across runs and all the endogenous intensities are shifted by a constant corresponding to reference intensities.
 - **FALSE** : No normalization is performed. If you had your own normalization before MSstats, you should use **Normalization=FALSE**.
 - NOTE : If there are multiple fractionations or injections for one sample, normalization is perform by each fractionation or different m/z range from multiple injections.
- **nameStandards** : Only for **Normalization='globalStandards'**, global standard peptide or Protein names, which you can assume that they have the same abundance across MS runs, should be assigned in the vector for this option.
- **featureSubset** :
 - ‘**all**’ : Use all features in the dataset.
 - ‘**top3**’ : Use top 3 features which have highest average of $\log_2(\text{intensity})$ across runs.
 - ‘**topN**’ : Use top N features which have highest average of $\log_2(\text{intensity})$ across runs. It needs the input for **n_top_feature** option (ex. **n_top_feature=5** for top 5 features).
 - ‘**highQuality**’ : Detect and flag uninformative features (as **Uninformative** in the **feature_quality** column) and outliers (as **TRUE** in the **is_outlier** column). These uninformative content may be excluded from run-level summarization by setting the **remove_uninformative_feature_outlier** option to **TRUE**.
- **summaryMethod** : Method for run-level summarization.
 - ‘**TMP**’ : Default. Tukey’s median polish (**medpolish** function in stats). Robust parameter estimation method with median across rows and columns.
 - ‘**linear**’ : Linear model (**lm** function). Average-based summarization.
- **MBimpute** : whether model-based imputation will be performed or not. Only for **summaryMethod='TMP'**.
 - **TRUE** : Default. Censored missing values will be imputed by Accelerated Failure Time model. Censored missing values will be determined by other options, **censoredInt** and **maxQuantileforCensored**
 - **FALSE** : No model-based imputation.
- **maxQuantileforCensored** : Maximum quantile for deciding censored missing value. Default is 0.999. If you don’t want to apply the threshold of noise intensity in your data, you can use **maxQuantileforCensored=NULL**.
- **censoredInt** : The processing tools report missing values differently. This option is for distinguish which value should be considered as missing, and further whether it is censored or at random.
 - ‘**NA**’ : Default. It assumes that all NAs in **Intensity** column are censored.
 - ‘**0**’ : It assumes that all values between 0 and 1 in **Intensity** column are censored. If there are NAs in **Intensity** with this option, NAs will be considered as random missing.

- **NULL** : It assumes that all missing values are randomly missing.
- **cutoffCensored** : cutoff value for AFT model. It is only with `censoredInt='NA'` or `censoredInt='0'`. If you have `censoredInt=NULL`, it assumes that there is no censored missing and any imputation will not be performed.
 - **'minFeature'** : cutoff for AFT model will be the minimum value for each feature across runs.
 - **'minRun'** : cutoff for AFT model will be the minimum value for each run across features.
 - **'minFeatureNRun'** : cutoff for AFT model will be the smallest value between minimum value of corresponding feature and minimum value of corresponding run.

A typical label-free DDA dataset may have many missing values and noisy features with outliers. **MSstats** supports several ways to deal with this. The default option for summarization is **TMP** (robust parameter estimation method with median across rows and columns) after imputation by **AFT** (accelerated failure time model, **MBimpute=TRUE**) based on censored intensity for NA (`censoredInt="NA"`) with a cutoff as the minimum value for a feature (`cutoffCensored="minFeature"`).

This process handles missing values through imputation and reduces the influence of the outliers using the **TMP** estimation. Note, however, that those runs with no measurements at all will be removed and not be used for any calculation.

default option

```
DDA2009.proposed <- dataProcess(raw = DDARawData,
                                normalization = 'equalizeMedians',
                                summaryMethod = 'TMP',
                                censoredInt = "NA",
                                cutoffCensored = "minFeature",
                                MBimpute = TRUE,
                                maxQuantileforCensored=0.999)

## ** Log2 intensities under cutoff = 13.456 were considered as censored missing values.
## ** Log2 intensities = NA were considered as censored missing values.
## ** Use all features that the dataset originally has.

##
## Summary of Features :
##
## # of Protein          count
## # of Peptides/Protein 11-32
## # of Transitions/Peptide 1-1
##
## Summary of Samples :
##
## C1 C2 C3 C4 C5 C6
## # of MS runs          3 3 3 3 3 3
## # of Biological Replicates 1 1 1 1 1 1
## # of Technical Replicates 3 3 3 3 3 3
##
## Summary of Missingness :
##
## # transitions are completely missing in at least one of the conditions : 90
##
## -> D.GPLTGTYR_23_23_NA_NA, F.HFWGSSDDQGSEHTVDR_402_402_NA_NA, G.PLTGTYR_8_8_NA_NA, H.SFNVEYDDSQ
##
## # run with 75% missing observations: 0
##
## == Start the summarization per subplot...
```

```
## |
##
## == the summarization per subplot is done.
```

Output of dataProcess Output of the dataProcess function contains the processed and run-level summarized data as well as relevant information for the summarization step.

```
# output of dataProcess includes several data types.
names(DDA2009.proposed)
```

```
## [1] "ProcessedData"      "RunlevelData"      "SummaryMethod"
## [4] "ModelQC"           "PredictBySurvival"
```

```
# the data after reformatting and normalization
head(DDA2009.proposed$ProcessedData)
```

```
##      PROTEIN      PEPTIDE TRANSITION
## 55   bovine      D.GPLTGTyr_23_23    NA_NA
## 937  bovine F.HFWGSSDDQGSEHTVDR_402_402 NA_NA
## 1628 bovine F.HWGSSDDQGSEHTVDR_229_229 NA_NA
## 19   bovine      G.PLTGTyr_8_8      NA_NA
## 1081 bovine      H.SFNVEYDDSQDK_465_465 NA_NA
## 469  bovine      K.AVVQDPALKPL_156_156 NA_NA
##
##      FEATURE LABEL GROUP_ORIGINAL SUBJECT_ORIGINAL
## 55      D.GPLTGTyr_23_23_NA_NA      L      C1      1
## 937 F.HFWGSSDDQGSEHTVDR_402_402_NA_NA L      C1      1
## 1628 F.HWGSSDDQGSEHTVDR_229_229_NA_NA L      C1      1
## 19      G.PLTGTyr_8_8_NA_NA      L      C1      1
## 1081 H.SFNVEYDDSQDK_465_465_NA_NA      L      C1      1
## 469 K.AVVQDPALKPL_156_156_NA_NA      L      C1      1
##
##      RUN GROUP SUBJECT INTENSITY SUBJECT_NESTED ABUNDANCE FRACTION originalRUN
## 55      1      1      1 757400.1      1.1 19.83052      1      1
## 937      1      1      1 2087125.8      1.1 21.29291      1      1
## 1628      1      1      1 1485145.8      1.1 20.80200      1      1
## 19      1      1      1 4986404.0      1.1 22.54939      1      1
## 1081      1      1      1 2488141.2      1.1 21.54646      1      1
## 469      1      1      1 7519322.0      1.1 23.14200      1      1
##
##      censored
## 55      FALSE
## 937      FALSE
## 1628      FALSE
## 19      FALSE
## 1081      FALSE
## 469      FALSE
```

DDA2009.TMP\$ProcessedData has the data after normalization and deciding the data-specific threshold for censored missing value. There are several new columns in the datasets. Also dataset is reformatted. **Intensity** column includes original intensities values in the input of dataProcess. **ABUNDANCE** column contains the *log2 transformed and normalized intensities and it will be used for run-level summarization*. **censored** column has the decision about censored missing or not, based on **censoredInt** and **maxQuantileforCensored** options. **ABUNDANCE** with TRUE value in **censored** column will be considered as censored missing and imputed with **MBimpute=TRUE** option. Censored missing will be distinguished in Profile plot from dataProcessPlots.

```
# run-level summarized data
head(DDA2009.proposed$RunlevelData)
```

```
##   RUN Protein LogIntensities NumMeasuredFeature MissingPercentage more50missing
## 1   1 bovine      21.28437             14      0.00000000      FALSE
## 2   2 bovine      20.85653             14      0.00000000      FALSE
## 3   3 bovine      20.67521             13      0.07142857      FALSE
## 4   4 bovine      21.60443             13      0.07142857      FALSE
## 5   5 bovine      21.82186             14      0.00000000      FALSE
## 6   6 bovine      21.20445             13      0.07142857      FALSE
##   NumImputedFeature originalRUN GROUP GROUP_ORIGINAL SUBJECT_ORIGINAL
## 1                0           1    1             C1             1
## 2                0           2    1             C1             1
## 3                1           3    1             C1             1
## 4                1           4    2             C2             1
## 5                0           5    2             C2             1
## 6                1           6    2             C2             1
##   SUBJECT_NESTED SUBJECT
## 1              1.1      1
## 2              1.1      1
## 3              1.1      1
## 4              2.1      1
## 5              2.1      1
## 6              2.1      1
```

DDA2009.TMP\$RunlevelData includes run-level summarized data based on DDA2009.TMP\$ProcessedData. LogIntensities is run-level summarized data and will be used for groupComparison function in next step. It will also be used for summarized profile plot (summaryPlot=TRUE for dataProcessPlots function with type='ProfilePlot'). NumMeasuredFeature shows how many features were used for summarization of the corresponding run and protein. MissingPercentage means the percentage of random and censored missing in the corresponding run and protein out of the total number of feature in the corresponding protein. more50missing means whether MissingPercentage is greater than 50% or not. NumImputedFeature shows how many features were imputed in the corresponding run and protein.


```
# here 'TMP' : It shows which summaryMethod is used for run-level summarization.
head(DDA2009.proposed$SummaryMethod)
```

```
## [1] "TMP"
```


4.1.3 Visualization of processed data

Quality control and normalization effects QC plot visualizes potential systematic biases between mass spectrometry runs. Also it can be used to assess the effects of the normalization step. After constant normalization, the median intensities of reference transitions across all proteins should be equal between runs. After quantile normalization, the distribution of reference intensities across all proteins should be equal between runs. This step generates two types of QC plots: one for all the proteins combined, and the other separately for each protein (produced in a separate pdf file). These plots can be generated for either all proteins at once or each protein individually if we have a large dataset. The example below shows both options.

```
# use type="QCplot" with all proteins
# change the upper limit of y-axis=35
# set up the size of pdf
dataProcessPlots(data = DDA2009.proposed, type="QCplot", ylimUp=35,
                  width=5, height=5)
```
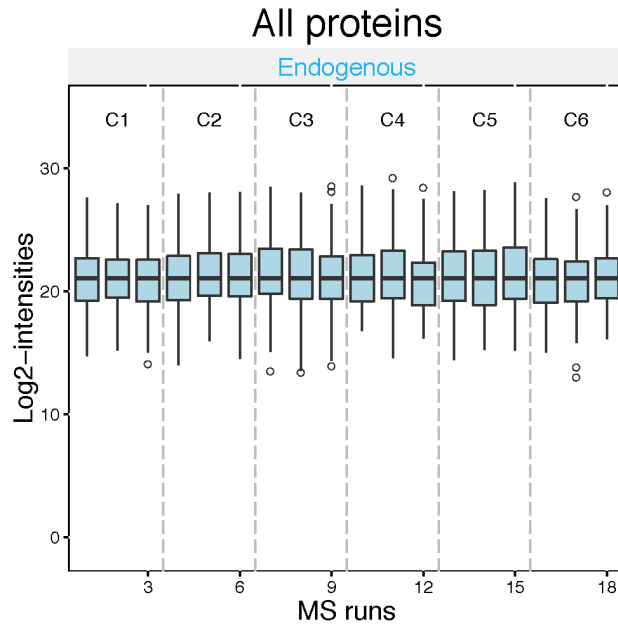
```
# use type="QCplot" for 1st protein only
# change the upper limit of y-axis=35
# set up the size of pdf
dataProcessPlots(data = DDA2009.proposed, type="QCplot", which.Protein=1,
  ylimUp=35, width=5, height=5)
```

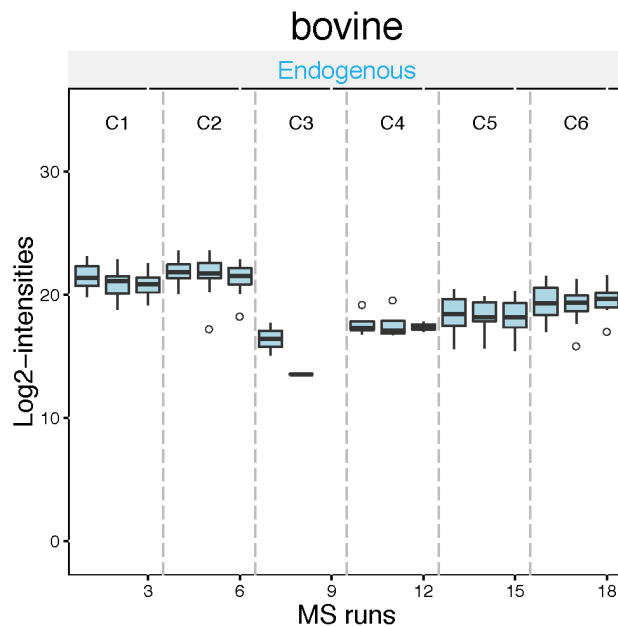


NOTE Don't worry about warning messages as below. It means NA values are not included in the plot, which is a proper way for this case.

Warning messages:
 1: Removed 698 rows containing non-finite values
 (stat_boxplot).
 ...


Profile plot Profile plot helps identify potential sources of variation (both variation of interest and nuisance variation) for each protein. Such plots should be done after the normalization. Profile plots with summarization present the effects of the summarization step by showing all individual measurements of a protein and their summarized intensity. With `type="profileplot"`, two pdfs will be generated. The first pdf includes plots (per protein) to show individual measurement for each peptide (peptide for DDA, transition for SRM or DIA) across runs, grouped per condition. **Each peptide has a different color/type layout.** Disconnected lines show that there are missing value (NA). To ignore these plots, please use the option `originalPlot=FALSE`. The second pdf, which is named with 'wSummarization' suffix, shows run-level summarized data per protein. The same peptides (or transition) in the first plot are presented in grey, with the summarized values (by TMP, in this example) overlaid in red. To ignore these plots with summarization, please use the option `summaryPlot=FALSE`.

```
dataProcessPlots(data = DDA2009.proposed, type="Profileplot", ylimUp=35,
                 featureName="NA", width=5, height=5, address="DDA2009_proposed_")
```

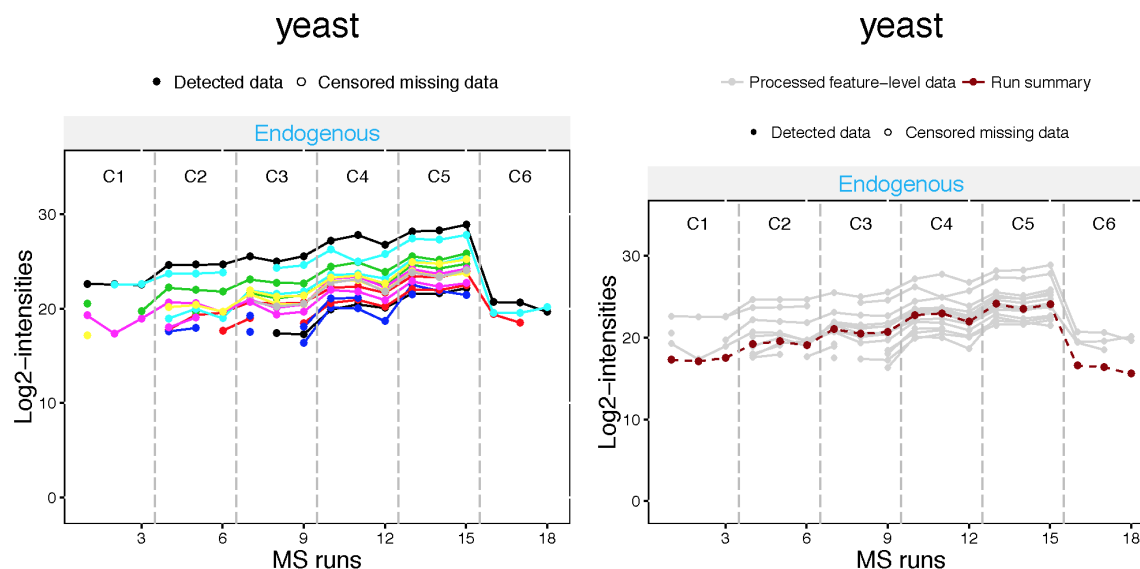

Condition plot Condition plot visualizes potential systematic differences in protein intensities between conditions. Dots indicate the mean of log2 intensities for each condition. With the option `interval='CI'` (default), error bars indicate the confidence interval with 0.95 significant level for each condition. With the option `interval='SD'`, error bars indicate the standard deviation among all feature intensities for each condition. **The intervals are for descriptive purposes only, as more refined model-based estimation is obtained as discussed below.** With the option `scale=TRUE`, the levels of conditions are scaled according to their labels. If `scale=FALSE` (default), the conditions on the x-axis are equally spaced.

```
dataProcessPlots(data = DDA2009.proposed, type="Conditionplot",
                 width=5, height=5, address="DDA2009_proposed_")
```

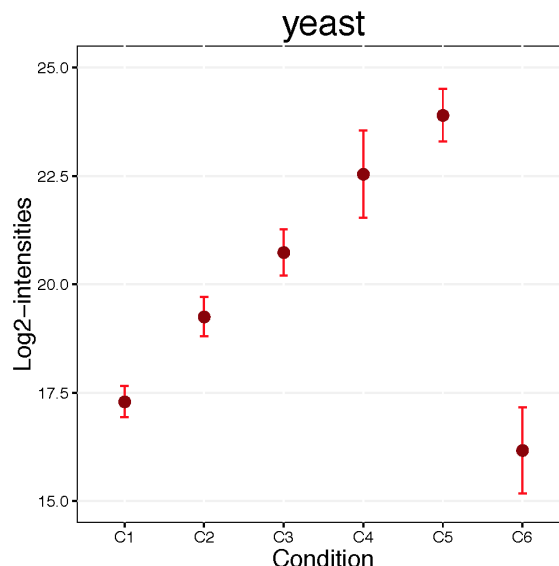

`dataProcessPlots` has a number of layout options, including size and description of axes labels, output file name etc for three types of plots above. The option `address` specifies the name of the folder storing pdf files with the plots. With the option `address=FALSE`, plots will be shown in the graphical window, but not saved in a file. If a file with this name already exists in working directory, a suffix with a number will be appended to the file name. In this way a record of all the analyses is kept.

For more details, visit the help file using the following code.

```
?dataProcessPlots
```

4.1.4 Different imputation options

Here is the summary of combinations for imputation options with `summaryMethod='TMP'`.

- `censoredInt=NULL` : It assumes that all intensities are missing at random and there is no action for missing values with `MBimpute=FALSE`. If you have `MBimpute=TRUE` with `censoredInt=NULL`, there will be error message to fix either `MBimpute` or `censoredInt` options.
- `censoredInt='NA' or '0' & MBimpute=TRUE` : AFT model-based imputation using `cutoffCensored` value in the AFT model.
- `censoredInt='NA' or '0' & MBimpute=FALSE` : censored intensities (here NA's) will be replaced with the value specified in `cutoffCensored`.

NOTE1 The default option for `cutoffCensored` is `minFeature`, taking the minimum value for the corresponding feature. With this option, those runs with substantial missing measurements will be biased by the cutoff value. In such case, you may remove the runs that have more than 50% missing values from the analysis with the option `remove50missing=TRUE`.

NOTE2 In case that there are completely missing measurements in a run for a protein, any imputation will not be performed. In addition, the condition, which has no measurement at all in a protein, will be not imputed.

Here is the example of `dataProcess` option without imputation, assuming that all missing values are random.

```
# No imputation
DDA2009.TMP <- dataProcess(raw = DDARawData,
                           normalization = 'equalizeMedians',
```

```

summaryMethod = 'TMP',
censoredInt = NULL, MBimpute=FALSE)

## ** Use all features that the dataset origianally has.
##
## Summary of Features :
## count
## # of Protein 6
## # of Peptides/Protein 11-32
## # of Transitions/Peptide 1-1
##
## Summary of Samples :
## C1 C2 C3 C4 C5 C6
## # of MS runs 3 3 3 3 3 3
## # of Biological Replicates 1 1 1 1 1 1
## # of Technical Replicates 3 3 3 3 3 3
##
## Summary of Missingness :
## # transitions are completely missing in at least one of the conditions : 90
## -> D.GPLTGYR_23_23_NA_NA, F.HFWGSSDDQGSEHTVDR_402_402_NA_NA, G.PLTGYR_8_8_NA_NA, H.SFNVEYDDSQ
##
## # run with 75% missing observations: 0
##
## == Start the summarization per subplot...
## |
##
## == the summarization per subplot is done.


```

These plots can be used compare and select among different options for imputation (e.g., TMP with or without considering missing values for summarization in dataProcess).

```

dataProcessPlots(data = DDA2009.TMP, type="Profileplot", ylimUp=35,
featureName="NA", width=5, height=5, address="DDA2009_TMP_")

```

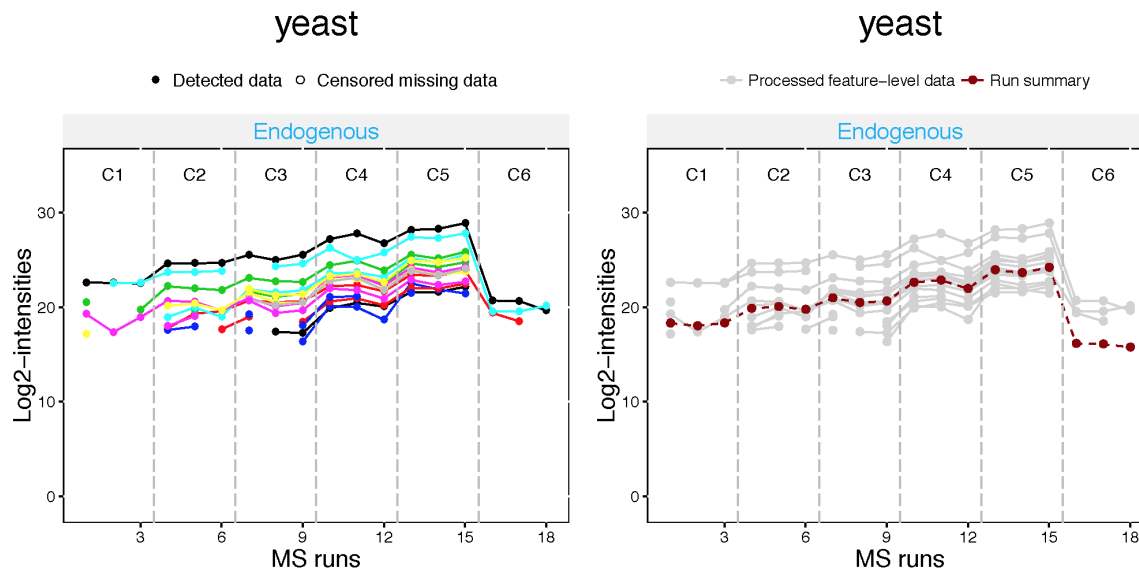

While original profile plots are the same, summarization plots reveal differences, especially for conditions 'C1' and 'C2' in 'yeast' protein, which have many missing values. Without imputation, summarized values in 'C1' group is higher than with imputation for missing values.

4.1.5 Feature selection

A feature selection module is integrated in the `dataProcess` function, and is performed with the option `featureSubset="highQuality"`. In the `ProcessedData` element of the output returned by `dataProcess`, the option adds two columns `feature_quality` and `is_outlier`, to highlight uninformative features that are inconsistent with the consensus profile ("Uninformative" in `feature_quality`) and outliers (TRUE in `is_outlier`). The uninformative features and outliers can be separately investigated, curated, or removed. The option `remove_uninformative_feature_outlier=TRUE` removes the detected uninformative features and outliers.

```
# Feature selection
DDA2009.inf <- dataProcess(raw = DDARawData,
                           normalization = 'equalizeMedians',
                           summaryMethod = 'TMP',
                           featureSubset = "highQuality",
                           remove_uninformative_feature_outlier = TRUE)
```

```
##
## Summary of Features :
## count
## # of Protein 6
## # of Peptides/Protein 11-32
## # of Transitions/Peptide 1-1
##
## Summary of Samples :
## C1 C2 C3 C4 C5 C6
## # of MS runs 3 3 3 3 3 3
## # of Biological Replicates 1 1 1 1 1 1
## # of Technical Replicates 3 3 3 3 3 3
## |
```

```
processed.inf <- DDA2009.inf$ProcessedData
```

```
table(processed.inf$feature_quality)
```

```
##
##    Informative Uninformative
##          2052           18
```

```
head(processed.inf[processed.inf$feature_quality == "Uninformative", ])
```

```
##      PROTEIN          PEPTIDE TRANSITION          FEATURE
## 87 cyc_horse A.PGFTYTDANKNK_367_367      NA_NA A.PGFTYTDANKNK_367_367_NA_NA
## 88 cyc_horse A.PGFTYTDANKNK_367_367      NA_NA A.PGFTYTDANKNK_367_367_NA_NA
## 89 cyc_horse A.PGFTYTDANKNK_367_367      NA_NA A.PGFTYTDANKNK_367_367_NA_NA
## 90 cyc_horse A.PGFTYTDANKNK_367_367      NA_NA A.PGFTYTDANKNK_367_367_NA_NA
## 91 cyc_horse A.PGFTYTDANKNK_367_367      NA_NA A.PGFTYTDANKNK_367_367_NA_NA
## 92 cyc_horse A.PGFTYTDANKNK_367_367      NA_NA A.PGFTYTDANKNK_367_367_NA_NA
##      LABEL GROUP_ORIGINAL SUBJECT_ORIGINAL RUN GROUP SUBJECT INTENSITY
## 87      L          C1              1  1  1      1 187598.06
## 88      L          C1              1  2  1      1      NA
## 89      L          C1              1  3  1      1  54184.12
## 90      L          C2              1  4  2      1 435188.97
## 91      L          C2              1  5  2      1 235653.19
## 92      L          C2              1  6  2      1      NA
##      SUBJECT_NESTED ABUNDANCE FRACTION originalRUN censored feature_quality
## 87              1.1  17.81711      1              1  FALSE  Uninformative
## 88              1.1      NA      1              2   TRUE  Uninformative
## 89              1.1  16.22112      1              3  FALSE  Uninformative
## 90              2.1  18.32737      1              4  FALSE  Uninformative
## 91              2.1  17.62688      1              5  FALSE  Uninformative
## 92              2.1      NA      1              6   TRUE  Uninformative
##      is_outlier
## 87      FALSE
## 88      FALSE
## 89      FALSE
## 90      FALSE
## 91      FALSE
## 92      FALSE
```


4.1.6 Different summarization options

Besides summarizing observations with the median polish method, MSstats also offers a summarization option using linear model with option `summaryMethod="linear"` with `censoredInt=NULL` assumes that all NA's are missing at random and uses `lm` for parameter estimation.

```
# linear model (lm) with run and feature
DDA2009.linear <- dataProcess(raw = DDARawData,
                             normalization = 'equalizeMedians',
                             summaryMethod = 'linear',
                             censoredInt = NULL,
                             MBimpute = FALSE)
```

Profile plots below can be used compare among different options for summarization (e.g., TMP with or without imputation vs `linear` for summarization in `dataProcess`).

```
dataProcessPlots(data = DDA2009.linear, type="Profileplot", ylimUp=35,
  featureName="NA", width=5, height=5, address="DDA2009_linear_")
```

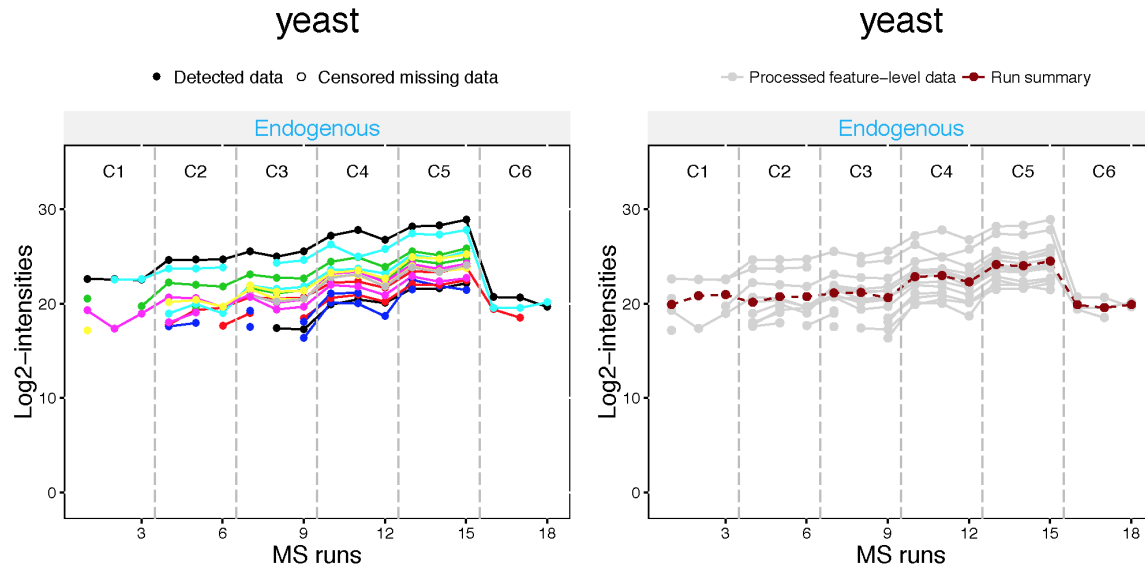

While original profile plots are the same, summarization plots reveal differences, especially for conditions 'C1', 'C2', and 'C6' in 'yeast' protein, which have many missing values. Summarized values with linear model in these groups are much higher than those with TMP considering missing values or not.

4.1.7 Finding differentially abundant proteins across conditions

Comparing conditions with groupComparison With the normalized data and run-level summarized data obtained by applying one of the `dataProcess` summarization methods, it is of general interest to find proteins changing between groups of conditions. Within `MSstats` this can be done by using the `groupComparison` function, which takes the output of the `dataProcess` function as input.

```
?groupComparison
```

In addition to the processed data, the `groupComparison` function requires a contrast matrix to define the comparison to be made. The contrast matrix is created with each condition in column and each comparison in row. Note that the conditions are arranged in **alphabetical order**. The order of condition that `MSstats` recognizes can be shown by using `levels`:

```
levels(DDA2009.TMP$ProcessedData$GROUP_ORIGINAL)
```

```
## [1] "C1" "C2" "C3" "C4" "C5" "C6"
```

Entries in each row of the contrast matrix are filled in with 0, 1, or -1 to specify the comparison, where 0 is for conditions we would like to ignore, 1 is for conditions we would like to put in the numerator of the ratio or fold-change, and -1 is for conditions we would like to put in the denominator of the ratio or fold-change.

For example, if you want to compare C2-C1, which means $\log(C2)-\log(C1)$ and the same as $\log(C2/C1)$, set '1' for C2 and '-1' for C1 in the row. Combining multiple groups for comparison is also possible. For example, if you want to compare between average of C2 and C3 and average of C1, $(C3+C2)/2-C1$ as formula, set '-1' for C1, '0.5' for C2 and '0.5' for C3, and '0' for rest of groups.

```
comparison1 <- matrix(c(-1,1,0,0,0,0),nrow=1)
comparison2 <- matrix(c(0,-1,1,0,0,0),nrow=1)
comparison3 <- matrix(c(0,0,-1,1,0,0),nrow=1)
```

```

comparison4 <- matrix(c(0,0,0,-1,1,0),nrow=1)
comparison5 <- matrix(c(0,0,0,0,-1,1),nrow=1)
comparison6 <- matrix(c(1,0,0,0,0,-1),nrow=1)

comparison<-rbind(comparison1,comparison2,comparison3,comparison4,comparison5,comparison6)
row.names(comparison) <- c("C2-C1","C3-C2","C4-C3","C5-C4","C6-C5","C1-C6")

```

With the contrast matrix specified, group comparison can be performed as follows.

```
DDA2009.comparisons <- groupComparison(contrast.matrix = comparison, data = DDA2009.proposed)
```

```
##      |
```

Output of the `groupComparison` function contains three data frames:

```

# output from groupComparison function has three data frames
names(DDA2009.comparisons)

```

```
## [1] "ComparisonResult" "ModelQC" "fittedmodel"
```

Results of the statistical comparison are stored in the data frame named `ComparisonResult`:

```

# name of columns in result data.frame
head(DDA2009.comparisons$ComparisonResult)

```

```

##      Protein Label      log2FC      SE      Tvalue DF      pvalue      adj.pvalue
## 1    bovine C2-C1  0.6048799 0.4245943  1.424607 11 1.820186e-01 1.820186e-01
## 2    chicken C2-C1  0.7876884 0.2205455  3.571545 12 3.841470e-03 4.609764e-03
## 3   cyc_horse C2-C1  1.1294964 0.1955787  5.775149 12 8.809149e-05 1.761830e-04
## 4  myg_horse C2-C1 -7.9717333 0.2807086 -28.398612 12 2.254641e-12 1.352785e-11
## 5     rabbit C2-C1  1.0617105 0.2209612  4.804963 12 4.298913e-04 6.448369e-04
## 6      yeast C2-C1  1.9575344 0.3134408  6.245309 12 4.284464e-05 1.285339e-04
##      issue MissingPercentage ImputationPercentage
## 1     NA          0.03571429          0.03571429
## 2     NA          0.25757576          0.25757576
## 3     NA          0.15104167          0.15104167
## 4     NA          0.45833333          0.45833333
## 5     NA          0.72043011          0.72043011
## 6     NA          0.56666667          0.56666667

```

The result of the test for differential abundance is a table with columns **Protein**, **Label** (of the comparison), **log2FC** (log2 fold change), **SE** (standard error of the log2 fold change), **Tvalue** (test statistic of the Student test), **DF** (degree of freedom of the Student test), **pvalue** (raw p-values), **adj.pvalue** (p-values adjusted among all the proteins in the specific comparison using the approach by Benjamini and Hochberg (Benjamini and Hochberg 1955)). The cutoff of the adjusted p-value corresponds to the cutoff of the False Discovery Rate (Benjamini and Hochberg 1955). The positive values of **log2FC** for **Label=C2-C1** indicate evidence in favor of $C2 > C1$ (i.e. proteins upregulated in C2), while the negative values indicate evidence in favor of $C2 < C1$ (i.e. proteins downregulated in C2), as compared to C1. The same model can be used to perform several comparisons of conditions simultaneously in the same protein.

NOTE `issue` column shows if there is any issue for inference in corresponding protein and comparison, for example, `OneConditionMissing` or `CompleteMissing`. If one of condition for comparison is completely missing, it would flag with `OneConditionMissing` with `adj.pvalue=0` and `log2FC=Inf` or `-Inf` even though `pvalue=NA`. For example, if you want to compare 'Condition A - Condition B', but condition B has complete missing, `log2FC=Inf` and `adj.pvalue=0`. `SE`, `Tvalue`, and `pvalue` will be `NA`. if you want to compare 'Condition A - Condition B', but condition A has complete missing, then `log2FC=-Inf` and `adj.pvalue=0`. But, please be careful for using this `log2FC` and `adj.pvalue`.

Based on the comparison results and desired significance level, a short list of the differentially abundant proteins can be obtained for further investigation:

```
# get only significant proteins and comparisons among all comparisons
# To simultaneously controll the overall FDR at the level, 0.05
SignificantProteins <- with(DDA2009.comparisons,
                           ComparisonResult[ComparisonResult$adj.pvalue < 0.05, ])
nrow(SignificantProteins)

## [1] 34
```

4.1.8 Verifying the assumption of the model


Results based on the statistical models are accurate as long as the assumptions of the models hold. Here we focus on the assumption of the Normal distribution of the measurement errors, and also on the assumption of constant variance of the measurement errors (if this option is specified in the model above). The assumptions can be checked by examining the residuals of the model fit (i.e., the deviations of the observed intensities of the transition from their model-based predictions).

`modelBasedQCPlots` function generates residual plots and Normal quantile-quantile plots for each protein, taking as input the results of model fitting and testing in `groupComparison`. Normal quantile-quantile plot with the option `type='QQPlots'` illustrates that such deviations from constant variance can be mistaken for deviations from Normality. Only large deviations of transition intensities from the straight line are problematic.

Residual plot with the option `type='ResidualPlots'` shows variance of the residuals that is associated with the mean feature intensity. Any specific pattern, such as increasing or decreasing by predicted abundance, is problematic.

```
# normal quantile-quantile plots
modelBasedQCPlots(data=DDA2009.comparisons, type="QQPlots",
                  width=5, height=5, address="DDA2009_proposed_")

# residual plots
modelBasedQCPlots(data=DDA2009.comparisons, type="ResidualPlots",
                  width=5, height=5, address="DDA2009_proposed_")
```

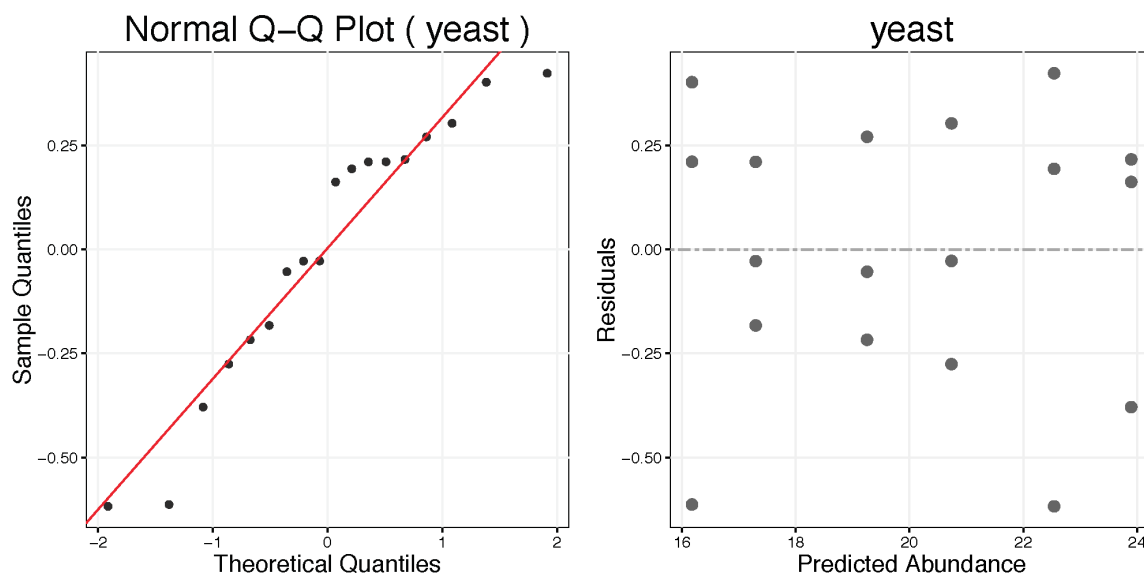

For more details, visit the help file using the following code.


```
?modelBasedQCPlots
```

4.1.9 Visualization of differentially abundant proteins

Volcano plots Volcano plots visualize the outcome of one comparison between conditions for all the proteins, and combine the information on statistical and practical significance. The y-axis displays the FDR-adjusted p-values on the negative log10 scale, representing statistical significance. The horizontal dashed line shows the FDR cutoff. The points above the FDR cutoff line are statistically significant proteins that are differentially abundant across conditions. These points are colored in red and blue for upregulated and downregulated proteins, respectively. The x-axis is the model-based estimate of fold change on log scale (the base of logarithm transform is the same as specified in the `logTrans` option of the `dataProcess` function), and represents practical significance. It is possible to specify a practical significance cutoff based on the estimate of fold change in addition to the statistical significance cutoff. If the fold change cutoff is specified, the points above the horizontal cutoff line but within the vertical cutoff line will be considered as not differentially abundant (and will be colored in black). The practical significance cutoff should only be applied in addition to the statistical significance cutoff (i.e., the fold change alone does not present enough evidence for differential abundance).

```
groupComparisonPlots(data = DDA2009.comparisons$ComparisonResult, type = 'VolcanoPlot',  
                     width=5, height=5, address="DDA2009_proposed_")
```

‘VolcanoPlot.pdf’ will be saved under the folder you assigned. It has the plots per comparison defined in `contrast.matrix`. Please check `?groupComparisonPlots` for detail, such as labelling protein names, size of dots, font sizes, etc. Below is one of volcano plots, for comparison ‘C1-C6’ including protein name labelling. Protein name will be shown for significant proteins, without overlapping protein names each other.

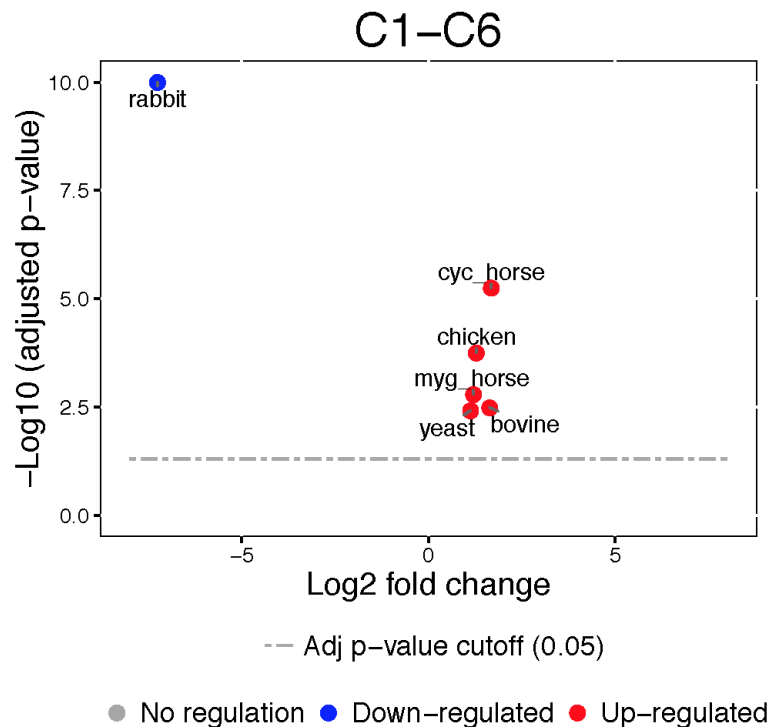


Heatmap Heatmaps illustrate the patterns of up- and down-regulation of proteins in several comparisons. Columns in the heatmaps are comparison of conditions assigned in `contrast.matrix`, and rows are proteins. The heatmaps display **signed FDR-adjusted p-values of the tests, colored in red/blue for significantly up-/down-regulated proteins, while taking into account the specified FDR cutoff and the additional optional fold change cutoff**. Brighter colors indicate stronger evidence in favor of differential abundance. Black color represents proteins that are not significantly differentially abundant.

NOTE To draw heatmap, at least two comparisons are needed.

The rows and columns of the heatmaps can be ordered with the option `clustering`, which performs hierarchical clustering with the Ward method (minimum variance). The option `clustering='protein'` (default) clusters the rows (proteins) in the space of comparisons, based on the values of $(\text{sign of comparison}) \cdot (-\log_2(\text{adjusted p-values}))$. The option `clustering='comparison'` clusters the columns in the space of proteins, based on the values of $(\text{sign of comparison}) \cdot (-\log_2(\text{adjusted p-value}))$. The option `clustering='both'` reorders both columns and rows.


```
groupComparisonPlots(data = DDA2009.comparisons$ComparisonResult, type = 'Heatmap')
```

‘Heatmap.pdf’ will be saved under the folder you assigned. Below is one example, showing the results for several comparisons simultaneously.

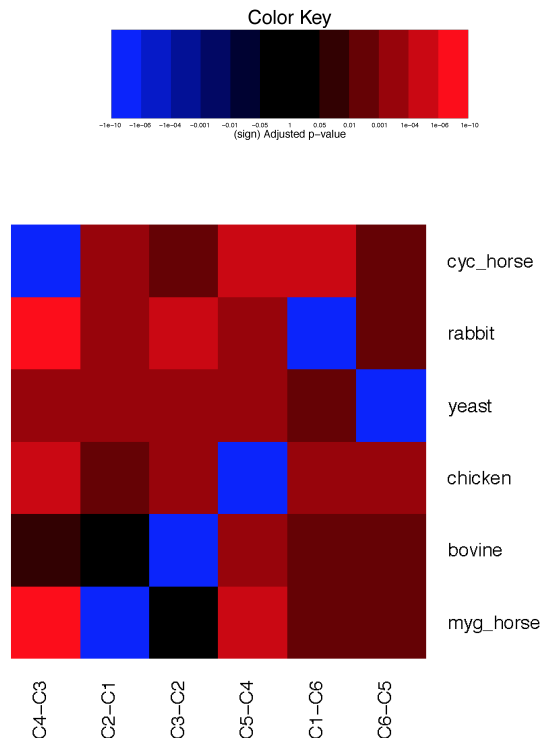

Comparison plots Comparison plots illustrate model-based estimates of log-fold changes, and the associated uncertainty, in several comparisons of conditions for one protein. X-axis is the comparison of interest. Y-axis is the log fold change. The dots are the model-based estimates of log-fold change, and the error bars are the model-based 95% confidence intervals (the option `sig` can be used to change the significance level of significance). For simplicity, the confidence intervals are adjusted for multiple comparisons within protein only, using the Bonferroni approach. For proteins with N comparisons, the individual confidence intervals are at the level of $1\text{-sig}/N$.

```
groupComparisonPlots(data=DDA2009.comparisons$ComparisonResult, type="ComparisonPlot",
                     width=5, height=5, address="DDA2009_proposed_")
```

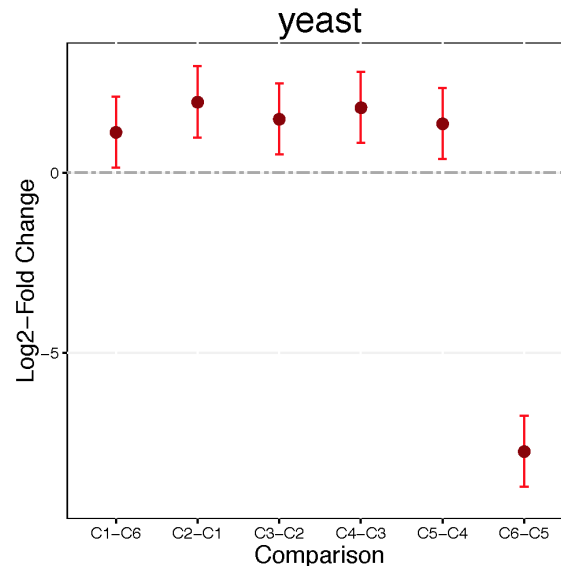

For further details, such as labelling protein names, size of dots, font sizes, etc., visit the help file using the following code.

```
?groupComparisonPlots
```

4.1.10 Sample size calculation for a future experiment

This last analysis step views the dataset as a pilot study of a future experiment, utilizes its variance components, and calculates the minimal number of replicates required in a future experiment to achieve the desired statistical power. The calculation is performed by the function `designSampleSize`, which takes as input the fitted model in `groupComparison`. Sample size calculation assumes same experimental design (i.e. group comparison, time course or paired design) as in the current dataset, and uses the model fit to estimate the median variance components across all the proteins. Finally, sample size calculation assumes that a large proportion of proteins (specifically, 99%) will not change in abundance in the future experiment. This assumption also provides conservative results. Using the estimated variance components, the function relates the number of biological replicates per condition (`numSample`, rounded to 0 decimal), average statistical power across all the proteins (`power`), minimal fold change that we would like to detect (can be specified as a range, e.g. `desiredFC=c(1.1, 2)`), and the False Discovery Rate (FDR). The user should specify all these quantities but one, and the function will solve for the remainder. The quantity to solve for should be set to `= TRUE`.

```
# Minimal number of biological replicates per condition
result.sample <- designSampleSize(data=DDA2009.comparisons$fittedmodel, numSample=TRUE,
                                   desiredFC=c(1.25, 3), FDR=0.05, power=0.8)
result.sample
```

##	desiredFC	numSample	FDR	power	CV
## 1	1.250	35	0.05	0.8	0.004
## 2	1.275	30	0.05	0.8	0.005
## 3	1.300	25	0.05	0.8	0.006
## 4	1.325	22	0.05	0.8	0.007
## 5	1.350	19	0.05	0.8	0.007
## 6	1.375	17	0.05	0.8	0.008

## 7	1.400	16 0.05	0.8 0.009
## 8	1.425	14 0.05	0.8 0.010
## 9	1.450	13 0.05	0.8 0.010
## 10	1.475	12 0.05	0.8 0.011
## 11	1.500	11 0.05	0.8 0.012
## 12	1.525	10 0.05	0.8 0.013
## 13	1.550	9 0.05	0.8 0.014
## 14	1.575	9 0.05	0.8 0.014
## 15	1.600	8 0.05	0.8 0.015
## 16	1.625	7 0.05	0.8 0.017
## 17	1.650	7 0.05	0.8 0.017
## 18	1.675	7 0.05	0.8 0.016
## 19	1.700	6 0.05	0.8 0.019
## 20	1.725	6 0.05	0.8 0.018
## 21	1.750	6 0.05	0.8 0.018
## 22	1.775	5 0.05	0.8 0.022
## 23	1.800	5 0.05	0.8 0.021
## 24	1.825	5 0.05	0.8 0.021
## 25	1.850	5 0.05	0.8 0.021
## 26	1.875	4 0.05	0.8 0.026
## 27	1.900	4 0.05	0.8 0.025
## 28	1.925	4 0.05	0.8 0.025
## 29	1.950	4 0.05	0.8 0.025
## 30	1.975	4 0.05	0.8 0.024
## 31	2.000	4 0.05	0.8 0.024
## 32	2.025	4 0.05	0.8 0.024
## 33	2.050	3 0.05	0.8 0.031
## 34	2.075	3 0.05	0.8 0.031
## 35	2.100	3 0.05	0.8 0.030
## 36	2.125	3 0.05	0.8 0.030
## 37	2.150	3 0.05	0.8 0.030
## 38	2.175	3 0.05	0.8 0.029
## 39	2.200	3 0.05	0.8 0.029
## 40	2.225	3 0.05	0.8 0.029
## 41	2.250	3 0.05	0.8 0.028
## 42	2.275	3 0.05	0.8 0.028
## 43	2.300	3 0.05	0.8 0.028
## 44	2.325	2 0.05	0.8 0.041
## 45	2.350	2 0.05	0.8 0.041
## 46	2.375	2 0.05	0.8 0.040
## 47	2.400	2 0.05	0.8 0.040
## 48	2.425	2 0.05	0.8 0.039
## 49	2.450	2 0.05	0.8 0.039
## 50	2.475	2 0.05	0.8 0.039
## 51	2.500	2 0.05	0.8 0.038
## 52	2.525	2 0.05	0.8 0.038
## 53	2.550	2 0.05	0.8 0.038
## 54	2.575	2 0.05	0.8 0.037
## 55	2.600	2 0.05	0.8 0.037
## 56	2.625	2 0.05	0.8 0.036
## 57	2.650	2 0.05	0.8 0.036
## 58	2.675	2 0.05	0.8 0.036
## 59	2.700	2 0.05	0.8 0.035
## 60	2.725	2 0.05	0.8 0.035

```
## 61      2.750      2 0.05    0.8 0.035
## 62      2.775      2 0.05    0.8 0.034
## 63      2.800      2 0.05    0.8 0.034
## 64      2.825      2 0.05    0.8 0.034
## 65      2.850      2 0.05    0.8 0.034
## 66      2.875      2 0.05    0.8 0.033
## 67      2.900      2 0.05    0.8 0.033
## 68      2.925      2 0.05    0.8 0.033
## 69      2.950      1 0.05    0.8 0.065
## 70      2.975      1 0.05    0.8 0.064
## 71      3.000      1 0.05    0.8 0.064
```

```
# Power calculation
```

```
result.power <- designSampleSize(data=DDA2009.comparisons$fittedmodel, numSample=3,
                                desiredFC=c(1.25, 3), FDR=0.05, power=TRUE)
result.power
```


```
##      desiredFC numSample  FDR power    CV
## 1      1.250      3 0.05  0.01 0.051
## 2      1.275      3 0.05  0.01 0.050
## 3      1.300      3 0.05  0.01 0.049
## 4      1.325      3 0.05  0.01 0.048
## 5      1.350      3 0.05  0.01 0.047
## 6      1.375      3 0.05  0.01 0.046
## 7      1.400      3 0.05  0.01 0.046
## 8      1.425      3 0.05  0.01 0.045
## 9      1.450      3 0.05  0.01 0.044
## 10     1.475      3 0.05  0.01 0.043
## 11     1.500      3 0.05  0.02 0.043
## 12     1.525      3 0.05  0.03 0.042
## 13     1.550      3 0.05  0.04 0.041
## 14     1.575      3 0.05  0.06 0.041
## 15     1.600      3 0.05  0.08 0.040
## 16     1.625      3 0.05  0.10 0.039
## 17     1.650      3 0.05  0.13 0.039
## 18     1.675      3 0.05  0.16 0.038
## 19     1.700      3 0.05  0.20 0.038
## 20     1.725      3 0.05  0.23 0.037
## 21     1.750      3 0.05  0.27 0.036
## 22     1.775      3 0.05  0.31 0.036
## 23     1.800      3 0.05  0.35 0.035
## 24     1.825      3 0.05  0.39 0.035
## 25     1.850      3 0.05  0.43 0.034
## 26     1.875      3 0.05  0.47 0.034
## 27     1.900      3 0.05  0.51 0.034
## 28     1.925      3 0.05  0.55 0.033
## 29     1.950      3 0.05  0.58 0.033
## 30     1.975      3 0.05  0.61 0.032
## 31     2.000      3 0.05  0.65 0.032
## 32     2.025      3 0.05  0.68 0.032
## 33     2.050      3 0.05  0.71 0.031
## 34     2.075      3 0.05  0.73 0.031
## 35     2.100      3 0.05  0.76 0.030
## 36     2.125      3 0.05  0.78 0.030
## 37     2.150      3 0.05  0.80 0.030
```

```
## 38      2.175      3 0.05  0.82 0.029
## 39      2.200      3 0.05  0.84 0.029
## 40      2.225      3 0.05  0.86 0.029
## 41      2.250      3 0.05  0.87 0.028
## 42      2.275      3 0.05  0.88 0.028
## 43      2.300      3 0.05  0.90 0.028
## 44      2.325      3 0.05  0.91 0.027
## 45      2.350      3 0.05  0.92 0.027
## 46      2.375      3 0.05  0.93 0.027
## 47      2.400      3 0.05  0.93 0.027
## 48      2.425      3 0.05  0.94 0.026
## 49      2.450      3 0.05  0.95 0.026
## 50      2.475      3 0.05  0.95 0.026
## 51      2.500      3 0.05  0.96 0.026
## 52      2.525      3 0.05  0.96 0.025
## 53      2.550      3 0.05  0.97 0.025
## 54      2.575      3 0.05  0.97 0.025
## 55      2.600      3 0.05  0.98 0.025
## 56      2.625      3 0.05  0.98 0.024
## 57      2.650      3 0.05  0.98 0.024
## 58      2.675      3 0.05  0.98 0.024
## 59      2.700      3 0.05  0.99 0.024
## 60      2.725      3 0.05  0.99 0.023
## 61      2.750      3 0.05  0.99 0.023
## 62      2.775      3 0.05  0.99 0.023
## 63      2.800      3 0.05  0.99 0.023
## 64      2.825      3 0.05  0.99 0.023
## 65      2.850      3 0.05  0.99 0.022
## 66      2.875      3 0.05  0.99 0.022
## 67      2.900      3 0.05  0.99 0.022
## 68      2.925      3 0.05  0.99 0.022
## 69      2.950      3 0.05  0.99 0.022
## 70      2.975      3 0.05  0.99 0.021
## 71      3.000      3 0.05  0.99 0.021
```

For further details, visit the help file using the following code.

```
?designSampleSize
```

Visualization of sample size calculations The calculated relationship between the number of biological replicates per condition (`numSample`), average statistical power across all the proteins (`power`), minimal fold change that we would like to detect (`desiredFC`), and the False Discovery Rate (`FDR`) can be visualized using the function `designSampleSizePlots`. The function takes as input the output of `designSampleSize`.

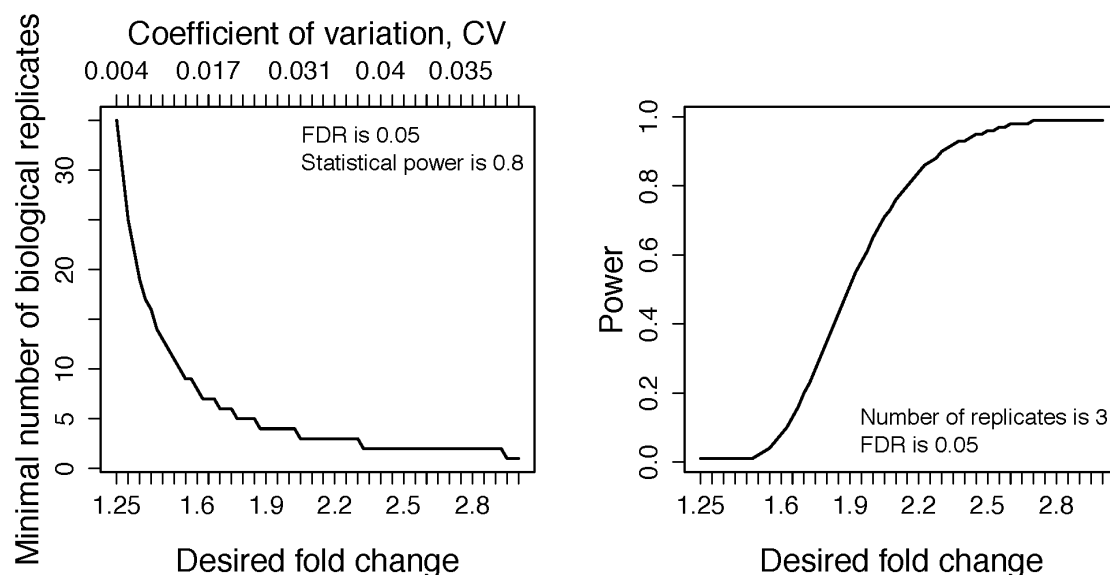

For further details, visit the help file using the following code.

```
?designSampleSizePlots
```

4.1.11 Quantification of protein abundance in individual samples or conditions

Many downstream analysis steps (such as clustering or classification of individual samples in the space of their protein profiles) require summary values of protein abundance in each biological replicate or in each condition, on a relative scale that is comparable between runs.

`dataProcess` function performs model-based run-level summarization. `quantification` function enables subject-level summarization or group-level summarization with the run-level summarization from `dataProcess`.

The option, `type='sample'` (default), performs sample quantification, i.e. it outputs the estimates of relative protein abundance in each biological replicate. If there are technical replicates for biological replicates, sample quantification will be the median among technical replicates. If there is no technical replicate for biological replicate (sample), sample quantification will be the same as run-level summarization. In presence of completely missing values in biological replicate, the estimates will be zero.

The option `type='group'` performs group quantification, i.e. it outputs the estimates of relative protein abundance in each condition, summarized over the biological replicates (median among sample quantification). In presence of completely missing values in a condition, the estimates will be zero.

MSstats supports two output formats. The option `format='matrix'` (default) outputs an array where rows are `proteins`, and columns are `conditions` (for group quantification), or combinations of biological replicate and condition ids (for sample quantification). The option `format='long'` produces an array where each row corresponding to relative protein abundances, and columns are `Protein`, `Condition`, `LogIntensities` (and `BioReplicate` in the case of sample quantification).

```
subQuant <- quantification(DDA2009.proposed)
head(subQuant)
```

```
##      Protein    C1_1    C2_1    C3_1    C4_1    C5_1    C6_1
## 1    bovine 20.85653 21.60443 14.32690 16.10441 17.63141 19.27802
## 2   chicken 18.48792 19.43204 20.41274 22.42284 15.92462 17.09803
## 3  cyc_horse 20.25927 21.33967 22.22028 15.85252 17.62720 18.45536
## 4  myg_horse 22.66495 14.73701 14.99667 18.61740 20.26392 21.52022
```



```
## 5    rabbit 14.89507 15.88492 17.43767 20.19014 21.27964 22.07550
## 6    yeast 17.26792 19.19987 20.71073 22.73666 24.06156 16.38660
```

```
groupQuant <- quantification(DDA2009.proposed, type='group')
head(groupQuant)
```

```
##      Protein      C1      C2      C3      C4      C5      C6
## 1    bovine 20.85653 21.60443 14.32690 16.10441 17.63141 19.27802
## 2   chicken 18.48792 19.43204 20.41274 22.42284 15.92462 17.09803
## 3 cyc_horse 20.25927 21.33967 22.22028 15.85252 17.62720 18.45536
## 4 myg_horse 22.66495 14.73701 14.99667 18.61740 20.26392 21.52022
## 5    rabbit 14.89507 15.88492 17.43767 20.19014 21.27964 22.07550
## 6    yeast 17.26792 19.19987 20.71073 22.73666 24.06156 16.38660
```

For further details, visit the help file using the following code.

```
?quantification
```

4.2 Suggested workflow with Skyline output for DDA

This section describes steps and considerations to properly format data processed by Skyline, prior to the MSstats analysis. In the following example, the raw files for the benchmark dataset (Choi, M. and Eren-Dogu, Z. F. and Colangelo, C. and Cottrell, J. and Hoopmann, M. R. and Kapp, E. A. and Kim, S. and Lam, H. and Neubert, T. A. and Palmblad, M. and Phinney, B. S. and Weintraub, S. T. and MacLean, B. and Vitek, O. 2017) are used. Dataset was processed and quantified with Skyline (3.5.0.9319). Details for data processing are described in Choi, et al., 2017 and Panorama Web <https://panoramaweb.org/iPRG-2015.url> for iProphet cut-off 0.15. The datasets and details for data processing are available in MassIVE.quant, MSV000079843, Reanalysis : RMSV000000249.1

4.2.1 Load Skyline output

This required input data is generated automatically when using MSstats report format in Skyline. We first load and access the dataset processed by Skyline. The name of saved file from Skyline using MSstats report format is 'Choi2017_DDA_Skyline_input.csv' under the folder named dda_skyline.

```
# Read output from skyline
raw <- read.csv("dda_skyline/Choi2017_DDA_Skyline_input.csv")
```

We can read csv file. Here we will load R data file which is the exactly same data in Choi2017_DDA_Skyline_input.csv file.

```
# Load R data, which is converted from csv file, output from skyline
load("dda_skyline/iprg.skyline.rda")
raw <- iprg.skyline
head(raw)
```

```
##      ProteinName PeptideSequence PeptideModifiedSequence
## 1 DECOY_sp|POCF18|YM085_YEAST    KDMYGNPFQK            KDM[+16] YGNPFQK
## 2 DECOY_sp|POCF18|YM085_YEAST    KDMYGNPFQK            KDM[+16] YGNPFQK
## 3 DECOY_sp|POCF18|YM085_YEAST    KDMYGNPFQK            KDM[+16] YGNPFQK
## 4 DECOY_sp|POCF18|YM085_YEAST    KDMYGNPFQK            KDM[+16] YGNPFQK
## 5 DECOY_sp|POCF18|YM085_YEAST    KDMYGNPFQK            KDM[+16] YGNPFQK
## 6 DECOY_sp|POCF18|YM085_YEAST    KDMYGNPFQK            KDM[+16] YGNPFQK
##      PrecursorCharge PrecursorMz FragmentIon ProductCharge ProductMz
```

```
## 1      3      415.1974 precursor      3 415.1974
## 2      3      415.1974 precursor      3 415.1974
## 3      3      415.1974 precursor      3 415.1974
## 4      3      415.1974 precursor      3 415.1974
## 5      3      415.1974 precursor      3 415.1974
## 6      3      415.1974 precursor      3 415.1974
##      IsotopeLabelType Condition BioReplicate      FileName
## 1      light Condition1      1 JD_06232014_sample1-A.raw
## 2      light Condition1      2 JD_06232014_sample1-B.raw
## 3      light Condition1      3 JD_06232014_sample1-C.raw
## 4      light Condition2      4 JD_06232014_sample2-A.raw
## 5      light Condition2      5 JD_06232014_sample2-B.raw
## 6      light Condition2      6 JD_06232014_sample2-C.raw
##      Area StandardType Truncated DetectionQValue
## 1 71765.046875      NA      False      #N/A
## 2 147327.265625      NA      False      #N/A
## 3  1373396.5      NA      False      #N/A
## 4 66387.4453125      NA      False      #N/A
## 5 107736.453125      NA      False      #N/A
## 6  380812.0625      NA      False      #N/A
```

Annotation information is required to fill in `Condition` and `BioReplicate` for corresponding `Run` information. Users have to prepare as csv or txt file like 'Choi2017_DDA_Skyline_annotation.csv', which includes `Run`, `Condition`, and `BioReplicate` information, and load it in R.

```
annot <- read.csv("dda_skyline/Choi2017_DDA_Skyline_annotation.csv", header=TRUE)
annot
```

```
##      Run Condition BioReplicate
## 1 JD_06232014_sample1-A.raw Condition1      1
## 2 JD_06232014_sample1-B.raw Condition1      1
## 3 JD_06232014_sample1-C.raw Condition1      1
## 4 JD_06232014_sample2-A.raw Condition2      2
## 5 JD_06232014_sample2-B.raw Condition2      2
## 6 JD_06232014_sample2-C.raw Condition2      2
## 7 JD_06232014_sample3-A.raw Condition3      3
## 8 JD_06232014_sample3-B.raw Condition3      3
## 9 JD_06232014_sample3-C.raw Condition3      3
## 10 JD_06232014_sample4-A.raw Condition4      4
## 11 JD_06232014_sample4-B.raw Condition4      4
## 12 JD_06232014_sample4-C.raw Condition4      4
```

4.2.2 Preprocessing with DDA experiment from Skyline output

The input data for `MSstats` is required to contain variables of `ProteinName`, `PeptideSequence`, `PrecursorCharge`, `FragmentIon`, `ProductCharge`, `IsotopeLabelType`, `Condition`, `BioReplicate`, `Run`, `Intensity`. These variable names should be fixed. `MSstats` input from Skyline adapts the column scheme of the dataset so that it fits `MSstats` input format. However there are several extra column names and also some of them need to be changed. `SkylinetoMSstatsFormat` function helps pre-processing for making right format of `MSstats` input from Skyline output. For example, it renames some column name, and replace truncated peak intensities with NA. Another important step is to handle isotopic peaks before using `dataProcess`. The output from Skyline for DDA experiment has several measurements of peak area from the monoisotopic, M+1 and M+2 peaks. To get a robust measure of peptide intensity, we can sum over isotopic peaks per peptide or use the highest peak. Here we take a summation per peptide ion.

Here is the summary of pre-processing steps in `SkylinetoMSstatsFormat` function.

SkylinetoMSstatsFormat

- Rename column names
- Remove decoy proteins and iRT peptides
- Remove shared peptides
- Replace NA for truncated intensities
- DDA: Sum of isotopic peaks per peptide and charge
- DIA : filter by q-value
- Remove features with all missing values or with less than 3 measurements across MS runs
- Remove protein with only one feature
- **Add annotation for experimental design** : Group, biological replicate, fraction information per MS run

Options for SkylinetoMSstatsFormat

- **annotation** : name of 'annotation.txt' or 'annotation.csv' data which includes Condition, BioReplicate, and Run. If annotation is already complete in Skyline, use `annotation=NULL` (default). It will use the annotation information from input.
- **removeiRT** : TRUE (default) will remove the proteins or peptides which are labeled 'iRT' in 'StandardType' column. FALSE will keep them.
- **filter_with_Qvalue** : TRUE (default) will filter out the intensities that have greater than `qvalue_cutoff` in `DetectionQValue` column. Those intensities will be replaced with **zero** and will be considered as censored missing values for imputation purpose.
- **qvalue_cutoff** : Cutoff for `DetectionQValue`. Default is 0.01.
- **useUniquePeptide** : TRUE (default) removes peptides that are assigned for more than one proteins. We assume to use unique peptide for each protein.
- **fewMeasurement** : remove or keep the feature with few measurements.
 - 'remove' : (default) remove the features that have 1 or 2 measurements across runs.
 - 'keep' : keep all the features. However, it could generate the error in the step for fitting the statistical model.
- **removeOxidationMpeptides** : TRUE will remove the peptides including 'oxidation (M)' in modification. FALSE is default.
- **removeProtein_with1Feature** : TRUE will remove the proteins which have only 1 peptide and charge. FALSE is default.

For further details, visit the help file using the following code.

```
?SkylinetoMSstatsFormat
```

Now, we use `SkylinetoMSstatsFormat` function for this example dataset. We chose to remove the proteins with only 1 peptide ion.

```
# reformatting and pre-processing for Skyline output.
quant <- SkylinetoMSstatsFormat(raw,
                                annotation = annot,
                                removeProtein_with1Feature = TRUE)
```

```
## ** Proteins, which names include DECOY, are removed.
```

```
## ** Peptides, that are used in more than one proteins, are removed.
## Warning in SkylinetoMSstatsFormat(raw, annotation = annot,
## removeProtein_with1Feature = TRUE): NAs introduced by coercion
## ** Truncated peaks are replaced with NA.
## ** For DDA datasets, three isotopic peaks per feature and run are summed.
## ** 4 features have all NAs or zero intensity values and are removed.
## ** 13 features have 1 or 2 intensities across runs and are removed.
## ** All proteins have at least two features.
```

This function shows the progress. The output of `SkylinetoMSstatsFormat`, called `quant`, is ready for next step.

```
head(quant)
```

```
##           ProteinName      PeptideSequence PrecursorCharge FragmentIon
## 1 sp|D6VTK4|STE2_YEAST EGEVEPVDM[+16] YTPDTAADEEARK          3         sum
## 2 sp|D6VTK4|STE2_YEAST EGEVEPVDM[+16] YTPDTAADEEARK          3         sum
## 3 sp|D6VTK4|STE2_YEAST EGEVEPVDM[+16] YTPDTAADEEARK          3         sum
## 4 sp|D6VTK4|STE2_YEAST EGEVEPVDM[+16] YTPDTAADEEARK          3         sum
## 5 sp|D6VTK4|STE2_YEAST EGEVEPVDM[+16] YTPDTAADEEARK          3         sum
## 6 sp|D6VTK4|STE2_YEAST EGEVEPVDM[+16] YTPDTAADEEARK          3         sum
##   ProductCharge IsotopeLabelType  Condition BioReplicate
## 1             NA                L Condition1             1
## 2             NA                L Condition2             2
## 3             NA                L Condition4             4
## 4             NA                L Condition2             2
## 5             NA                L Condition4             4
## 6             NA                L Condition3             3
##           Run Intensity StandardType
## 1 JD_06232014_sample1_C.raw 7863713.2      NA
## 2 JD_06232014_sample2_A.raw 977615.1      NA
## 3 JD_06232014_sample4_B.raw 4102785.2      NA
## 4 JD_06232014_sample2_C.raw 6547298.9      NA
## 5 JD_06232014_sample4_C.raw 3972463.8      NA
## 6 JD_06232014_sample3_B.raw 8896050.8      NA
```

4.2.3 Different options for Skyline in `dataProcess`

The difference between output from Skyline and other spectral processing tool is that Skyline distinguishes random missing (NA) by technical issues and low noisy intensity due to less than limit of detection. The output from Skyline can have both NA (expect small number of NAs or none of them) and very small intensity close to zero (less than 1 in intensity) and those should be treated different types of missing. In `dataProcess`, users need to use `censoredInt='0'` for Skyline output, which means to distinguish between NA as random missing and 0 as censored missing.

```
skyline.proposed <- dataProcess(quant,
                                normalization='equalizeMedian',
                                summaryMethod="TMP",
                                cutoffCensored="minFeature",
                                censoredInt="0", ## !! important
                                MBimpute=TRUE,
                                maxQuantileforCensored=0.999)
```

Further steps is the same as in general workflow (section 4.1).

4.3 Suggested workflow with MaxQuant output for DDA

The following R code chunks show steps to format a MaxQuant output for analysis by **MSstats**. In the following example, the raw files for the benchmark dataset (Choi, M. and Eren-Dogu, Z. F. and Colangelo, C. and Cottrell, J. and Hoopmann, M. R. and Kapp, E. A. and Kim, S. and Lam, H. and Neubert, T. A. and Palmblad, M. and Phinney, B. S. and Weintraub, S. T. and MacLean, B. and Vitek, O. 2017) are used. MS/MS spectra were searched using MaxQuant (v1.5.1.2) and Andromeda as search engine. The datasets and details for data processing are available in MassIVE.quant, MSV000079843, Reanalysis : RMSV000000249.2

4.3.1 Load MaxQuant outputs

Three files should be prepared before **MSstats**. Two files, ‘proteinGroups.txt’ and ‘evidence.txt’ are outputs from MaxQuant.

```
## First, get protein ID information
proteinGroups <- read.table("dda_maxquant/Choi2017_DDA_MaxQuant_proteinGroups.txt", sep = "\t", header = TRUE)

## Read in MaxQuant file: evidence.txt
infile <- read.table("dda_maxquant/Choi2017_DDA_MaxQuant_evidence.txt", sep = "\t", header = TRUE)
```

One file is for annotation information, required to fill in **Condition** and **BioReplicate** for corresponding Run information. Users have to prepare as csv or txt file like ‘Choi2017_DDA_MaxQuant_annotation.csv’, which includes Run, Condition, and BioReplicate information, and load it in R.

```
## Read in annotation including condition and biological replicates: annotation.csv
annot <- read.csv("dda_maxquant/Choi2017_DDA_MaxQuant_annotation.csv", header = TRUE)
annot
```

##		Raw.file	Condition	BioReplicate	Experiment	IsotopeLabelType
## 1	JD_06232014_sample1-A	Condition1	1	sample1_A	L	
## 2	JD_06232014_sample2_A	Condition2	2	sample2_A	L	
## 3	JD_06232014_sample4_B	Condition4	4	sample4_B	L	
## 4	JD_06232014_sample1_B	Condition1	1	sample1_B	L	
## 5	JD_06232014_sample1_C	Condition1	1	sample1_C	L	
## 6	JD_06232014_sample2_B	Condition2	2	sample2_B	L	
## 7	JD_06232014_sample2_C	Condition2	2	sample2_C	L	
## 8	JD_06232014_sample3_A	Condition3	3	sample3_A	L	
## 9	JD_06232014_sample3_B	Condition3	3	sample3_B	L	
## 10	JD_06232014_sample3_C	Condition3	3	sample3_C	L	
## 11	JD_06232014_sample4-A	Condition4	4	sample4_A	L	
## 12	JD_06232014_sample4_C	Condition4	4	sample4_C	L	

4.3.2 Preprocessing with DDA experiment from MaxQuant output

MaxQtoMSstatsFormat function helps pre-processing for making right format of **MSstats** input from MaxQuant output. Basically, this function gets peptide ion intensity from ‘evidence.txt’ file. In addition, there are several steps to filter out or to modify the data in order to get required information.

Here is the summary of pre-processing steps in **MaxQtoMSstatsFormat** function.

MaxQtoMSstatsFormat

- Remove '+' contaminant, reverse, Only.identified.by.site proteins
- Use 'protein.IDs' in proteinGroups.txt
- Extract essential information(columns)
- Remove shared peptides
- Aggregate multiple measurement per feature and run
- Remove features with all missing values or with less than 3 measurements across MS runs
- Remove protein with only one feature
- **Add annotation for experimental design** : Group, biological replicate, fraction information per MS run

Options for MaxQtoMSstatsFormat

- **evidence** : name of 'evidence.txt' data, which includes feature-level data
- **proteinGroups** : name of 'proteinGroups.txt' data. It needs to matching protein group ID. If proteinGroups=NULL, use 'Proteins' column in 'evidence.txt'.
- **annotation** :name of 'annotation.txt' or 'annotation.csv' data which includes Raw.file, Condition, BioReplicate, Run, and IsotopeLabelType information.
- **proteinID** : which column in evidence.txt will be used for ProteinName in MSstats.
 - 'Proteins' : (default) Proteins column will be used.
 - 'Leading.razor.protein' : Leading.razor.protein column will be used.
- **useUniquePeptide** : TRUE (default) removes peptides that are assigned for more than one proteins. We assume to use unique peptide for each protein.
- **summaryforMultipleRows** : max(default), sum, or mean. MSstats assumes that there is only one measurement (peak intensity) for one feature and one run. When there are multiple measurements for certain feature and certain run, MSstats need to know which measurements need to be used for further analysis. Users can use highest(max), sum or mean among multiple measurements for one feature and one run.
- **fewMeasurement** : remove or keep the featurew with few measurements.
 - 'remove' : (default) remove the features that have 1 or 2 measurements across runs.
 - 'keep' : keep all the features. However, it could generate the error in the step for fitting the statistical model.
- **removeMpeptides** : TRUE will remove the peptides including 'M' sequence. FALSE is default.
- **removeOxidationMpeptides** : TRUE will remove the peptides including 'oxidation (M)' in modification. FALSE is default.
- **removeProtein_with1Peptide** : TRUE will remove the proteins which have only 1 peptide and charge. FALSE is default.

For further details, visit the help file using the following code.

```
## check options for converting format
?MaxQtoMSstatsFormat
```

Now, we use MaxQtoMSstatsFormat function for this example dataset. We chose to remove the proteins with only 1 peptide ion.

```
quant <- MaxQtoMSstatsFormat(evidence=infile, annotation=annot, proteinGroups=proteinGroups,
                             removeProtein_with1Peptide=TRUE)
```

** + Contaminant, + Reverse, + Only.identified.by.site, proteins are removed.
 ## ** Peptide and charge, that have 1 or 2 measurements across runs, are removed.
 ## ** 282 proteins, which have only peptide and charge in a protein, are removed among 3157 proteins.
 This function shows the progress. The output of MaxQtoMSstatsFormat, called quant, is ready for next step.

```
## now 'quant' is ready for MSstats
head(quant)
```

##	ProteinName	PeptideSequence	PrecursorCharge	FragmentIon	ProductCharge
## 1	D6VTK4	EGEVEPVDMYTPDTAADEEARK	3	NA	NA
## 2	D6VTK4	FYPGTLSSFQTD SINNDK	2	NA	NA
## 3	D6VTK4	IGPFADASYK	2	NA	NA
## 4	D6VTK4	NQFYQLPTPTSSK	2	NA	NA
## 5	D6VTK4	TFVSETADDIEK	2	NA	NA
## 6	D6VTK4	TNTITSDFTTSTDR	2	NA	NA

##	IsotopeLabelType	Condition	BioReplicate	Run	Intensity
## 1	L	Condition1	1	JD_06232014_sample1_B	87141000
## 2	L	Condition1	1	JD_06232014_sample1_B	46167000
## 3	L	Condition1	1	JD_06232014_sample1_B	45425000
## 4	L	Condition1	1	JD_06232014_sample1_B	47094000
## 5	L	Condition1	1	JD_06232014_sample1_B	NA
## 6	L	Condition1	1	JD_06232014_sample1_B	62786000

4.3.3 Different options for MaxQuant in dataProcess

MaxQuant has certain or fixed threshold for intensity value internally as an parameter. Intensities less than the threshold are reported as NA. All missing values are NA in output from MaxQuant. In dataProcess, users need to use censoredInt='NA'. Users can use the same choice for other options.

```
maxquant.proposed <- dataProcess(quant,
                                 normalization='equalizeMedian',
                                 summaryMethod="TMP",
                                 cutoffCensored="minFeature",
                                 censoredInt="NA", ## !! important for MaxQuant
                                 MBimpute=TRUE,
                                 maxQuantileforCensored=0.999)
```

Further steps is the same as in general workflow (section 4.1).

4.4 Suggested workflow with Progenesis output for DDA

This section describes steps and considerations to properly format data processed by Progenesis, prior to the MSstats analysis. In the following example, the raw files for the benchmark dataset (Choi, M. and Eren-Dogu, Z. F. and Colangelo, C. and Cottrell, J. and Hoopmann, M. R. and Kapp, E. A. and Kim, S. and Lam, H. and Neubert, T. A. and Palmblad, M. and Phinney, B. S. and Weintraub, S. T. and MacLean, B. and Vitek, O. 2017) are used. Peptide features were identified with the Progenesis algorithm (v4.0.6403), aligned across all files, and annotated with the peptide identification resulting from the database search result

from Comet. The datasets and details for data processing are available in MassIVE.quant, MSV000079843, Reanalysis : RMSV000000249.3

4.4.1 Load Progenesis output

Here is the expected input for MSstats, which is output of Progenesis.

```
## First, read output of Progenesis
```

```
raw <- read.csv("dda_progenesis/Choi2017_DDA_Progenesis_input.csv")
```

```
head(raw)
```

```
##      X      X.1    X.2      X.3      X.4
## 1
## 2      # Retention time (min) Charge      m/z      Measured mass
## 3      16      52.5563333333333      2 501.781277638303 1001.54800234285
## 4      32      38.15255      2 474.251481407549 946.488409881339
## 5 11167      36.2224333333333      2 474.25154745893 946.488541984099
## 6      41      45.5598      2 371.731536419815 741.448519905869
##      X.5      X.6    X.7      X.8      X.9
## 1
## 2      Mass error (u) Mass error (ppm) Score Sequence Modifications
## 3 -0.00255665715405939 -2.55269904358308      1 TANDVLITR
## 4 -0.00219111866147159 -2.31499251990099      1 VTDGVMVAR
## 5 -0.0020590159006133 -2.17542139186367      1 VTDGVMVAR
## 6 -0.00120309413080122 -1.62262402086192 0.9996 AGLNIVR
##      X.10
## 1
## 2      Accession
## 3 sp|P00549|KPYK1_YEAST
## 4 sp|P00549|KPYK1_YEAST
## 5 sp|P00549|KPYK1_YEAST
## 6 sp|P00549|KPYK1_YEAST
##
## 1
## 2
## 3 Pyruvate kinase 1 OS=Saccharomyces cerevisiae (strain ATCC 204508 \ S288c) GN=CDC19 PE=1 SV=2
## 4 Pyruvate kinase 1 OS=Saccharomyces cerevisiae (strain ATCC 204508 \ S288c) GN=CDC19 PE=1 SV=2
## 5 Pyruvate kinase 1 OS=Saccharomyces cerevisiae (strain ATCC 204508 \ S288c) GN=CDC19 PE=1 SV=2
## 6 Pyruvate kinase 1 OS=Saccharomyces cerevisiae (strain ATCC 204508 \ S288c) GN=CDC19 PE=1 SV=2
##      X.12      X.13      X.14
## 1
## 2 Use in quantitation Max fold change Highest mean condition
## 3      False 1.23101575731737      A
## 4      False 1.35108253622201      B
## 5      False 1.25419527606242      B
## 6      False 1.04868912680216      A
##      X.15      X.16      X.17
## 1
## 2 Lowest mean condition Anova      Maximum CV
## 3      C 0.0522715538027003 12.0175133289667
## 4      A 0.0393818452091522 26.8776079679151
## 5      A 0.253277920596793 27.2310093101224
## 6      B 0.993981434364646 27.0631636013386
##      Normalized.abundance      X.18      X.19
```



```

## 1 A
## 2 JD_06232014_sample1-A JD_06232014_sample2_A JD_06232014_sample3_A
## 3 234646642.659118 246323351.490501 306102714.66799
## 4 179120293.733639 104309665.701784 136741892.392964
## 5 2233197.90782367 1134566.5998162 1574437.81004362
## 6 123797188.716029 122761256.64621 116107425.243685
## X.20 X.21 X.22
## 1 B
## 2 JD_06232014_sample4-A JD_06232014_sample1_B JD_06232014_sample2_B
## 3 257629531.217182 235779468.539422 236753257.546934
## 4 105011188.469111 182469696.644175 183285243.781685
## 5 1701362.72342001 2336796.51015815 1788630.29256942
## 6 63610598.879437 108255803.660911 108785457.069653
## X.23 X.24 X.25
## 1 C
## 2 JD_06232014_sample3_B JD_06232014_sample4_B JD_06232014_sample1_C
## 3 186699807.218591 242959514.796972 223435557.783206
## 4 162853464.030243 180957229.609825 186073547.948691
## 5 1932732.32106691 2274168.76697097 2182403.4051037
## 6 92469286.5619254 96974519.2450418 115737794.475905
## X.26 X.27 X.28
## 1
## 2 JD_06232014_sample2_C JD_06232014_sample3_C JD_06232014_sample4_C
## 3 220456628.684641 197285091.954229 207473304.500931
## 4 163946644.909844 150247529.397207 176820166.306319
## 5 2040841.8048229 1497501.50492545 2156318.09540588
## 6 100428622.768589 109442962.863466 83998024.0208531
## Raw.abundance X.29 X.30
## 1 A
## 2 JD_06232014_sample1-A JD_06232014_sample2_A JD_06232014_sample3_A
## 3 244531299.931508 221199440.186087 277078923.760572
## 4 186665863.932395 93670533.1411598 123776414.130536
## 5 2327272.96336292 1018845.73758267 1425154.0840074
## 6 129012233.635811 110240142.001841 105098448.610706
## X.31 X.32 X.33
## 1 B
## 2 JD_06232014_sample4-A JD_06232014_sample1_B JD_06232014_sample2_B
## 3 213112377.670857 265826760.55748 265610928.042007
## 4 86865756.2312952 205723291.596601 205625739.64841
## 5 1407375.36397932 2634593.46237989 2006645.04833363
## 6 52618990.9526949 122051719.672593 122045231.85501
## X.34 X.35 X.36
## 1 C
## 2 JD_06232014_sample3_B JD_06232014_sample4_B JD_06232014_sample1_C
## 3 219812880.452275 242959514.796972 210648874.118526
## 4 191737150.420336 180957229.609825 175425002.929311
## 5 2275521.67817463 2274168.76697097 2057509.66730008
## 6 108869636.960823 96974519.2450418 109114396.746851
## X.37 X.38 X.39
## 1
## 2 JD_06232014_sample2_C JD_06232014_sample3_C JD_06232014_sample4_C
## 3 207260386.701525 217212520.118758 194785326.277837
## 4 154133015.755399 165423774.187439 166006773.109084
## 5 1918679.71576989 1648761.5589421 2024448.99975338

```

```
## 6      94417098.3431607      120497608.498225      78861145.7987456
##      Spectral.counts      X.40      X.41
## 1      A
## 2 JD_06232014_sample1-A JD_06232014_sample2_A JD_06232014_sample3_A
## 3      2      1      1
## 4      1      1      1
## 5      1      0      1
## 6      2      2      2
##      X.42      X.43      X.44
## 1      B
## 2 JD_06232014_sample4-A JD_06232014_sample1_B JD_06232014_sample2_B
## 3      3      2      3
## 4      1      1      1
## 5      1      0      0
## 6      2      2      2
##      X.45      X.46      X.47
## 1      C
## 2 JD_06232014_sample3_B JD_06232014_sample4_B JD_06232014_sample1_C
## 3      2      2      2
## 4      1      1      1
## 5      0      1      1
## 6      2      2      2
##      X.48      X.49      X.50
## 1
## 2 JD_06232014_sample2_C JD_06232014_sample3_C JD_06232014_sample4_C
## 3      3      1      2
## 4      1      1      1
## 5      0      1      0
## 6      2      2      2
```

One file is for annotation information, required to fill in `Condition` and `BioReplicate` for corresponding Run information. Users have to prepare as csv or txt file like ‘Choi2017_DDA_Progenesis_annotation.csv’, which includes Run, Condition, and BioReplicate information, and load it in R.

```
## Read in annotation including condition and biological replicates
annot <- read.csv("dda_progenesis/Choi2017_DDA_Progenesis_annotation.csv", header = TRUE)
annot
```

```
##      Run Condition BioReplicate
## 1 JD_06232014_sample1-A Condition1      1
## 2 JD_06232014_sample2_A Condition2      2
## 3 JD_06232014_sample4_B Condition4      4
## 4 JD_06232014_sample1_B Condition1      1
## 5 JD_06232014_sample1_C Condition1      1
## 6 JD_06232014_sample2_B Condition2      2
## 7 JD_06232014_sample2_C Condition2      2
## 8 JD_06232014_sample3_A Condition3      3
## 9 JD_06232014_sample3_B Condition3      3
## 10 JD_06232014_sample3_C Condition3      3
## 11 JD_06232014_sample4-A Condition4      4
## 12 JD_06232014_sample4_C Condition4      4
```

4.4.2 Preprocessing with DDA experiment from progenesis output

The output from Progenesis includes peptide ion-level quantification for each MS runs. `ProgenesitoMSstatsFormat` function helps pre-processing for making right format of MSstats input from Progenesis output. Basically, this function reformats wide format to long format. It provide 'Raw.abundance', 'Normalized.abundance' and 'Spectral count' columns. This converter uses 'Raw.abundance' columns for Intensity values. In addition, there are several steps to filter out or to modify the data in order to get required information.

Here is the summary of pre-processing steps in `ProgenesitoMSstatsFormat` function.

ProgenesitoMSstatsFormat

- Extract essential information(columns)
- Subset the rows with `use in quantitation = true`
- Rename column names
- Remove empty protein and peptide sequence
- Remove duplicated measurements
- Remove shared peptides
- Aggregate multiple measurements per feature and run
- Remove features with all missing values or with less than 3 measurements across MS runs
- Remove protein with only one feature
- **Add annotation for experimental design** : Group, biological replicate, fraction information per MS run

Options for ProgenesitoMSstatsFormat

- **input** : name of Progenesis output, which is wide-format. 'Accession', 'Sequence', 'Modification', 'Charge' and one column for each run are required. 'Accession' column is used for ProteinName. 'Raw.abundance' is used for Intensity.
- **annotation** :name of 'annotation.txt' or 'annotation.csv' data which includes Condition, BioReplicate, and Run information. It will be matched with the column name of input for MS runs.
- **useUniquePeptide** : TRUE (default) removes peptides that are assigned for more than one proteins. We assume to use unique peptide for each protein.
- **summaryforMultipleRows** : max(default), sum, or mean. MSstats assumes that there is only one measurement (peak intensity) for one feature and one run. When there are multiple measurements for certain feature and certain run, MSstats need to know which measurements need to be used for further analysis. Users can use highest(max), sum or mean among multiple measurements for one feature and one run.
- **fewMeasurement** : remove or keep the featurew with few measurements.
 - 'remove' : (default) remove the features that have 1 or 2 measurements across runs.
 - 'keep' : keep all the features. However, it could generate the error in the step for fitting the statistical model.
- **removeOxidationMpeptides** : TRUE will remove the peptides including 'oxidation (M)' in modification. FALSE is default.
- **removeProtein_with1Peptide** : TRUE will remove the proteins which have only 1 peptide and charge. FALSE is default.

For further details, visit the help file using the following code.

```
## check options for converting format
?ProgenesisToMSstatsFormat
```

Now, we use `ProgenesisToMSstatsFormat` function for this example dataset. We chose to remove the proteins only 1 peptide ion.

```
quant <- ProgenesisToMSstatsFormat(raw, annotation=annot,
                                   removeProtein_with1Peptide = TRUE)
```

This function shows the progress. The output of `ProgenesisToMSstatsFormat`, called `quant`, is ready for next step.

```
## now 'quant' is ready for MSstats
head(quant)
```

##	ProteinName	PeptideSequence	PrecursorCharge	FragmentIon	ProductCharge
## 1	D6VTK4	EGEVEPVDMYTPDTAADEEARK	3	NA	NA
## 2	D6VTK4	FYPGTLSSFQTD SINNDK	2	NA	NA
## 3	D6VTK4	IGPFADASYK	2	NA	NA
## 4	D6VTK4	NQFYQLPTPTSSK	2	NA	NA
## 5	D6VTK4	TFVSETADDIEK	2	NA	NA
## 6	D6VTK4	TNTITSDFTTSTDR	2	NA	NA

##	IsotopeLabelType	Condition	BioReplicate	Run	Intensity
## 1	L	Condition1	1	JD_06232014_sample1_B	87141000
## 2	L	Condition1	1	JD_06232014_sample1_B	46167000
## 3	L	Condition1	1	JD_06232014_sample1_B	45425000
## 4	L	Condition1	1	JD_06232014_sample1_B	47094000
## 5	L	Condition1	1	JD_06232014_sample1_B	NA
## 6	L	Condition1	1	JD_06232014_sample1_B	62786000

4.4.3 Different options for Progenesis in dataProcess

Progenesis reports 0(zero) for missing values and does not have NA. Therefore, in `dataProcess`, users need to use `censoredInt='0'`. Users can use the same choice for other options.

```
progenesis.proposed <- dataProcess(quant,
                                   normalization='equalizeMedian',
                                   summaryMethod="TMP",
                                   cutoffCensored="minFeature",
                                   censoredInt="0", ## !! important
                                   MBimpute=TRUE,
                                   maxQuantileforCensored=0.999)
```

Further steps is the same as in general workflow (section 4.1).

4.5 Suggested workflow with Proteome Discoverer output for DDA

This section describes steps and considerations to properly format data processed by Proteome Discoverer, prior to the `MSstats` analysis. In the following example, another spike-in dataset processed by Proteome Discoverer is used to demonstrate. The datasets and details for data processing are available in `MassIVE.quant`, MSV000084181, Reanalysis : RMSV000000261.4

4.5.1 Load Proteome Discoverer output

The output from Proteome Discoverer includes several level of datasets. PSM sheet should be saved as csv as below. Here is the expected input for MSstats.

```
## Read PSM-level data
```

```
raw <- read.csv("dda_PD/ControlMixture_DDA_ProteomeDiscoverer_input.csv")
```

```
head(raw)
```

```
## Confidence.Level Search.ID Processing.Node.No Sequence Unique.Sequence.ID
## 1 High A 4 AALGVLR 2
## 2 High A 4 NLLLVK 4
## 3 High A 4 LIVVEK 5
## 4 High A 4 LLVDLK 6
## 5 High A 4 IITLLK 9
## 6 High A 4 HEFLR 10
## PSM.Ambiguity
## 1 Unambiguous
## 2 Unambiguous
## 3 Unambiguous
## 4 Unambiguous
## 5 Unambiguous
## 6 Unambiguous
## Protein.Descrip
## 1 Glycine--tRNA ligase beta subunit OS=Escherichia coli (strain K12) GN=glyS PE=1 SV=4 - [SYGB_E
## 2 50S ribosomal protein L3 OS=Escherichia coli (strain K12) GN=rplC PE=1 SV=1 - [RL3_E
## 3 50S ribosomal protein L4 OS=Escherichia coli (strain K12) GN=rplD PE=1 SV=1 - [RL4_E
## 4 Peptidyl-prolyl cis-trans isomerase D OS=Escherichia coli (strain K12) GN=ppiD PE=1 SV=1 - [PPID_E
## 5 3-dehydroquinate synthase OS=Escherichia coli (strain K12) GN=aroB PE=1 SV=1 - [AROB_E
## 6 GTP cyclohydrolase 1 OS=Escherichia coli (strain K12) GN=folE PE=1 SV=2 - [GCH1_E
## X..Proteins X..Protein.Groups Protein.Group.Accessions Modifications
## 1 1 1 P00961
## 2 1 1 P60438
## 3 1 1 P60723
## 4 1 1 POADY1
## 5 1 1 P07639
## 6 1 1 POA6T5
## Activation.Type DeltaScore DeltaCn Rank Search.Engine.Rank Precursor.Area
## 1 CID 1.0000 0 1 1 3.77e+07
## 2 CID 0.5455 0 1 1 6.59e+08
## 3 CID 0.0000 0 1 1 3.83e+08
## 4 CID 0.4062 0 1 1 1.42e+07
## 5 CID 1.0000 0 1 1 3.93e+07
## 6 CID 1.0000 0 1 1 2.80e+07
## QuanResultID Decoy.Peptides.Matched Exp.Value Homology.Threshold
## 1 NA 11 0.00033 13
## 2 NA 6 0.00940 13
## 3 NA 17 0.20000 13
## 4 NA 4 0.01300 13
## 5 NA NA 0.00860 13
## 6 NA 7 0.27000 13
## Identity.High Identity.Middle IonScore Peptides.Matched X..Missed.Cleavages
## 1 13 13 48 5 0
## 2 13 13 33 11 0
## 3 13 13 20 19 0
```

```

## 4      13      13      32      6      0
## 5      13      13      34      5      0
## 6      13      13      19      4      0
## Isolation.Interference.... Ion.Inject.Time..ms. Intensity Charge m.z..Da.
## 1      53      4 1700000      2 350.2295
## 2      8      2 2520000      2 350.2417
## 3     38      5  739000      2 350.7340
## 4     34      3 1520000      2 350.7342
## 5     13      2 2480000      2 350.7520
## 6     41     70  53500      2 351.1900
## MH...Da. Delta.Mass..Da. Delta.Mass..PPM. RT..min. First.Scan Last.Scan
## 1 699.4517      0      0.68 32.17      8180      8180
## 2 699.4761      0     -0.44 38.77     10907     10907
## 3 700.4607      0      0.41 27.49      6221      6221
## 4 700.4611      0      0.93 43.27     12766     12766
## 5 700.4968      0     -0.03 42.75     12552     12552
## 6 701.3728      0     -0.25 17.39      2693      2693
## MS.Order Ions.Matched Matched.Ions Total.Ions
## 1      MS2      Jun-50      6      50
## 2      MS2      May-52      5      52
## 3      MS2      May-40      5      40
## 4      MS2      May-40      5      40
## 5      MS2      Apr-40      4      40
## 6      MS2      Apr-32      4      32
## Spectrum.File Annotation
## 1 121219_S_CCES_01_01_LysC_Try_1to10_Mixt_1_1.raw      NA
## 2 121219_S_CCES_01_01_LysC_Try_1to10_Mixt_1_1.raw      NA
## 3 121219_S_CCES_01_01_LysC_Try_1to10_Mixt_1_1.raw      NA
## 4 121219_S_CCES_01_01_LysC_Try_1to10_Mixt_1_1.raw      NA
## 5 121219_S_CCES_01_01_LysC_Try_1to10_Mixt_1_1.raw      NA
## 6 121219_S_CCES_01_01_LysC_Try_1to10_Mixt_1_1.raw      NA

```

One file is for annotation information, required to fill in `Condition` and `BioReplicate` for corresponding `Run` information. Users have to prepare as csv or txt file like 'ControlMixture_DDA_ProteomeDiscoverer_annotation.csv', which includes `Run`, `Condition`, and `BioReplicate` information, and load it in R.

```

## Read in annotation including condition and biological replicates
annot <- read.csv("dda_PD/ControlMixture_DDA_ProteomeDiscoverer_annotation.csv", header = TRUE)
annot

```

```

## Run Condition BioReplicate
## 1 121219_S_CCES_01_01_LysC_Try_1to10_Mixt_1_1.raw Condition1 1
## 2 121219_S_CCES_01_02_LysC_Try_1to10_Mixt_1_2.raw Condition1 1
## 3 121219_S_CCES_01_03_LysC_Try_1to10_Mixt_1_3.raw Condition1 1
## 4 121219_S_CCES_01_04_LysC_Try_1to10_Mixt_2_1.raw Condition2 2
## 5 121219_S_CCES_01_05_LysC_Try_1to10_Mixt_2_2.raw Condition2 2
## 6 121219_S_CCES_01_06_LysC_Try_1to10_Mixt_2_3.raw Condition2 2
## 7 121219_S_CCES_01_07_LysC_Try_1to10_Mixt_3_1.raw Condition3 3
## 8 121219_S_CCES_01_08_LysC_Try_1to10_Mixt_3_2.raw Condition3 3
## 9 121219_S_CCES_01_09_LysC_Try_1to10_Mixt_3_3.raw Condition3 3
## 10 121219_S_CCES_01_10_LysC_Try_1to10_Mixt_4_1.raw Condition4 4
## 11 121219_S_CCES_01_11_LysC_Try_1to10_Mixt_4_2.raw Condition4 4
## 12 121219_S_CCES_01_12_LysC_Try_1to10_Mixt_4_3.raw Condition4 4
## 13 121219_S_CCES_01_13_LysC_Try_1to10_Mixt_5_1.raw Condition5 5
## 14 121219_S_CCES_01_14_LysC_Try_1to10_Mixt_5_2.raw Condition5 5

```

4.5.2 Preprocessing with DDA experiment from Proteome Discoverer output

PDtoMSstatsFormat function helps pre-processing for making right format of MSstats input from Proteome Discoverer output. `Protein.Group.Accessions` is used for `ProteinName`. The combination of `Sequence` and `Modifications` is used for `PeptideSequence`. `Charge` is used for `PrecursorCharge`. `Precursor.Area` is used for `Intensity`. In addition, there are several steps to filter out or to modify the data in order to get required information.

Here is the summary of pre-processing steps in PDtoMSstatsFormat function.

PDtoMSstatsFormat

- Extract essential information(columns)
- Rename column names
- Remove shared peptides
- Aggregate multiple measurements per feature and run : max or mean
- Remove features with all missing values or with less than 3 measurements across MS runs
- Remove protein with only one feature(Peptide ion)
- **Add annotation for experimental design** : Group, biological replicate, fraction information per MS run

Options for PDtoMSstatsFormat

- **input** : name of Proteome discover PSM output, which is long-format. “Protein.Group.Accessions”, “#Proteins”, “Sequence”, “Modifications”, “Charge”, “Intensity”, “Spectrum.File” are required.
- **annotation** : name of ‘annotation.txt’ or ‘annotation.csv’ data which includes Condition, BioReplicate, and Run information. ‘Run’ will be matched with ‘Spectrum.File’.
- **useNumProteinsColumn** : TRUE removes peptides which have more than 1 in # Proteins column of PD output.
- **useUniquePeptide** : TRUE (default) removes peptides that are assigned for more than one proteins. We assume to use unique peptide for each protein.
- **summaryforMultipleRows** : max(default), sum, or mean. MSstats assumes that there is only one measurement (peak intensity) for one feature and one run. When there are multiple measurements for certain feature and certain run, MSstats need to know which measurements need to be used for further analysis. Users can use highest(max), sum or mean among multiple measurements for one feature and one run.
- **fewMeasurement** : remove or keep the feature with few measurements.
 - ‘remove’ : (default) remove the features that have 1 or 2 measurements across runs.
 - ‘keep’ : keep all the features. However, it could generate the error in the step for fitting the statistical model.
- **removeOxidationMpeptides** : TRUE will remove the modified peptides including ‘Oxidation (M)’ in ‘Modifications’ column. FALSE is default.
- **removeProtein_with1Peptide** : TRUE will remove the proteins which have only 1 peptide and charge. FALSE is default.

- **which.quantification** : Use 'Precursor.Area'(default) column for quantified intensities. 'Intensity' or 'Area' can be used instead.
- **which.proteinid** : Use 'Protein.Accessions'(default) column for protein name. 'Master.Protein.Accessions' can be used instead.
- **which.sequence** : Use 'Sequence'(default) column for peptide sequence. 'Annotated.Sequence' can be used instead.

For further details, visit the help file using the following code.

```
## check options for converting format
?PDtoMSstatsFormat
```

Now, we use PDtoMSstatsFormat function for this example dataset. We chose to remove the proteins with only 1 peptide ion.

```
quant <- PDtoMSstatsFormat(raw,
                           annotation=annot,
                           which.proteinid = 'Protein.Group.Accessions',
                           removeProtein_with1Peptide=TRUE)
```

```
## ** Multiple measurements in a feature and a run are summarized by summaryforMultipleRows.
## ** 633 features have all NAs or zero intensity values and are removed.
## ** 241 proteins, which have only one feature in a protein, are removed among 1512 proteins.
```

This function shows the progress. The output of PDtoMSstatsFormat, called quant, is ready for next step.

```
head(quant)
```

```
## ProteinName PeptideModifiedSequence PrecursorCharge FragmentIon ProductCharge
## 1 P00961 AALGVLR_ 2 NA NA
## 2 P60438 NLLLVK_ 2 NA NA
## 3 P60723 LIVVEK_ 2 NA NA
## 4 POADY1 LLVDLK_ 2 NA NA
## 5 P07639 IITLLK_ 2 NA NA
## 6 POA6T5 HEFLR_ 2 NA NA
## IsotopeLabelType Condition BioReplicate
## 1 L Condition1 1
## 2 L Condition1 1
## 3 L Condition1 1
## 4 L Condition1 1
## 5 L Condition1 1
## 6 L Condition1 1
## Run Intensity
## 1 121219_S_CCES_01_01_LysC_Try_1to10_Mixt_1_1.raw 3.77e+07
## 2 121219_S_CCES_01_01_LysC_Try_1to10_Mixt_1_1.raw 6.59e+08
## 3 121219_S_CCES_01_01_LysC_Try_1to10_Mixt_1_1.raw 3.83e+08
## 4 121219_S_CCES_01_01_LysC_Try_1to10_Mixt_1_1.raw 1.42e+07
## 5 121219_S_CCES_01_01_LysC_Try_1to10_Mixt_1_1.raw 3.93e+07
## 6 121219_S_CCES_01_01_LysC_Try_1to10_Mixt_1_1.raw 2.80e+07
```

4.5.3 Different options for Proteome Discoverer in dataProcess

Progenesis reports NA for missing values. Therefore, in dataProcess, users need to use censoredInt='NA'. Users can use the same choice for other options.


```
pd.proposed <- dataProcess(quant,
  normalization='equalizeMedian',
  summaryMethod="TMP",
  cutoffCensored="minFeature",
  censoredInt="NA", ## !! important
  MBimpute=TRUE,
  maxQuantileforCensored=0.999)
```

Further steps is the same as in general workflow (section 4.1).

4.6 Suggested workflow with OpenMS output for DDA

This section describes steps and considerations to properly format data processed by OpenMS, prior to the MSstats analysis. In the following example, the raw files for the benchmark dataset (Choi, M. and Eren-Dogu, Z. F. and Colangelo, C. and Cottrell, J. and Hoopmann, M. R. and Kapp, E. A. and Kim, S. and Lam, H. and Neubert, T. A. and Palmblad, M. and Phinney, B. S. and Weintraub, S. T. and MacLean, B. and Vitek, O. 2017) are used. The datasets and details for data processing are available in OpenMS webpage, iPRG 2015 example dataset

4.6.1 Load OpenMS output

Here is the expected input for MSstats, which is output of OpenMS.

```
## Read PSM-level data
raw <- read.csv("dda_openms/iPRG_lfq_input_MSstats.csv")
head(raw)
```

##	ProteinName	PeptideSequence	PrecursorCharge	FragmentIon
## 1	sp P09938 RIR2_YEAST	AAADALSDLEIK	2	NA
## 2	sp P09938 RIR2_YEAST	AAADALSDLEIK	2	NA
## 3	sp P09938 RIR2_YEAST	AAADALSDLEIK	2	NA
## 4	sp P09938 RIR2_YEAST	AAADALSDLEIK	2	NA
## 5	sp P09938 RIR2_YEAST	AAADALSDLEIK	2	NA
## 6	sp P09938 RIR2_YEAST	AAADALSDLEIK	2	NA

##	ProductCharge	IsotopeLabelType	Condition	BioReplicate	Run	Intensity
## 1	0	L	1	1	1	391797000
## 2	0	L	4	10	10	103656000
## 3	0	L	4	11	11	361107000
## 4	0	L	1	2	2	456756000
## 5	0	L	1	3	3	389268000
## 6	0	L	2	4	4	433488000

If you follow the workflow in OpenMS KNIME, annotation information should be already filled in.

4.6.2 Preprocessing with DDA experiment from OpenMS output

OpenMStoMSstatsFormat function helps pre-processing for making right format of MSstats input from OpenMS output.

Here is the summary of pre-processing steps in OpenMStoMSstatsFormat function.

OpenMStoMSstatsFormat

- Extract essential information(columns)
- Rename column names
- Remove shared peptides
- Aggregate multiple measurements per feature and run
- Remove features with all missing values or with less than 3 measurements across MS runs
- Remove protein with only one feature
- **Add annotation for experimental design** : Group, biological replicate, fraction information per MS run

Options for OpenMStoMSstatsFormat

- **input** : name of MSstats input report from OpenMS, which includes feature(peptide ion)-level data.
- **annotation** : name of 'annotation.txt' data which includes Raw.file, Condition, BioReplicate, and Run information. If annotation is already complete in OpenMS, use annotation=NULL (default). It will use the annotation information from input.
- **useUniquePeptide** : TRUE (default) removes peptides that are assigned for more than one proteins. We assume to use unique peptide for each protein.
- **summaryforMultipleRows** : max(default), sum, or mean. MSstats assumes that there is only one measurement (peak intensity) for one feature and one run. When there are multiple measurements for certain feature and certain run, MSstats need to know which measurements need to be used for further analysis. Users can use highest(max), sum or mean among multiple measurements for one feature and one run.
- **fewMeasurement** : remove or keep the feature with few measurements.
 - **'remove'** : (default) remove the features that have 1 or 2 measurements across runs.
 - **'keep'** : keep all the features. However, it could generate the error in the step for fitting the statistical model.
- **removeProtein_with1Feature** : TRUE will remove the proteins which have only 1 peptide and charge. FALSE is default.

For further details, visit the help file using the following code.

```
## check options for converting format
?OpenMStoMSstatsFormat
```

Now, we use `OpenMStoMSstatsFormat` function for this example dataset. We chose to remove the proteins with shared peptides or only 1 peptide ion.

```
quant <- OpenMStoMSstatsFormat(raw,
                               removeProtein_with1Feature = TRUE)

## ** 0 features have all NAs or zero intensity values and are removed.
## ** All peptides are unique peptides in proteins.
## ** 909 features have 1 or 2 intensities across runs and are removed.
## ** 698 proteins, which have only one feature in a protein, are removed among 2538 proteins.
## ** No multiple measurements in a feature and a run.
```

```
## Warning in OpenMSstoMSstatsFormat(raw, removeProtein_with1Feature = TRUE): NAs
## introduced by coercion
```

This function shows the progress. The output of `OpenMSstoMSstatsFormat`, called `quant`, is ready for next step.

```
## now 'quant' is ready for MSstats
head(quant)
```

```
##           ProteinName      PeptideSequence PrecursorCharge FragmentIon
## 1 sp|D6VTK4|STE2_YEAST EGEVEPVDMYTPDTAADEEARK           3         NA
## 2 sp|D6VTK4|STE2_YEAST  FYPGTLSSSFQTD SINND AK           2         NA
## 3 sp|D6VTK4|STE2_YEAST          IGP FADASYK           2         NA
## 4 sp|D6VTK4|STE2_YEAST      NQFYQLPTPTSSK           2         NA
## 5 sp|D6VTK4|STE2_YEAST      TFVSETADDIEK           2         NA
## 6 sp|D6VTK4|STE2_YEAST      TNTITSDFTTSTDR           2         NA
##   ProductCharge IsotopeLabelType Condition BioReplicate Run Intensity
## 1             0                L          1             1  1  64757900
## 2             0                L          1             1  1  38852700
## 3             0                L          1             1  1  73225800
## 4             0                L          1             1  1  63139900
## 5             0                L          1             1  1         NA
## 6             0                L          1             1  1  58905300
```

4.6.3 Different options for OpenMS in dataProcess

Progenesis reports NA for missing values. Therefore, in `dataProcess`, users need to use `censoredInt='NA'`. Users can use the same choice for other options.

```
openms.proposed <- dataProcess(quant,
                               normalization='equalizeMedian',
                               summaryMethod="TMP",
                               cutoffCensored="minFeature",
                               censoredInt="NA",
                               MBimpute=TRUE,
                               maxQuantileforCensored=0.999)
```

Further steps is the same as in general workflow (section 4.1).

5. DIA analysis with MSstats

5.1 Suggested workflow with Skyline output for DIA

The analysis for DIA with Skyline output is the same as the workflow with Skyline output for DDA. Please check section 4.2 and options from `SkylinetoMSstatsFormat`.

5.2 Suggested workflow with Spectronaut output for DIA

This section describes steps and considerations to properly format data processed by Spectronaut for SWATH/DIA experiments, prior to the `MSstats` analysis. In the following example, the raw files for the

benchmark dataset (Bruderer et al. 2015) are quantified by Spectronaut. The datasets and details for data processing are available in MassIVE.quant, MSV000081828, Reanalysis : RMSV000000252.2

5.2.1 Load Spectronaut output

We first load and access the dataset processed by Spectronaut.

```
# Read output from Spectronaut
raw <- read.csv("dia_spectronaut/Bruderer2015_DIA_Spectronaut_input.xls", sep="\t")
```

One file is for annotation information, required to fill in `Condition` and `BioReplicate` for corresponding Run information. Users have to prepare as csv or txt file like 'Bruderer2015_DIA_Spectronaut_annotation.csv', which includes Run, Condition, and BioReplicate information, and load it in R.

```
## Read in annotation including condition and biological replicates
annot <- read.csv("dia_spectronaut/Bruderer2015_DIA_Spectronaut_annotation.csv", header = TRUE)
annot
```

##		Run	Condition	BioReplicate
## 1	B_D140314_SGSDSsample1_R01_MHRM		S1	S1
## 2	B_D140314_SGSDSsample1_R02_MHRM		S1	S1
## 3	B_D140314_SGSDSsample1_R03_MHRM		S1	S1
## 4	B_D140314_SGSDSsample2_R01_MHRM		S2	S2
## 5	B_D140314_SGSDSsample2_R02_MHRM		S2	S2
## 6	B_D140314_SGSDSsample2_R03_MHRM		S2	S2
## 7	B_D140314_SGSDSsample3_R01_MHRM		S3	S3
## 8	B_D140314_SGSDSsample3_R02_MHRM		S3	S3
## 9	B_D140314_SGSDSsample3_R03_MHRM		S3	S3
## 10	B_D140314_SGSDSsample4_R01_MHRM		S4	S4
## 11	B_D140314_SGSDSsample4_R02_MHRM		S4	S4
## 12	B_D140314_SGSDSsample4_R03_MHRM		S4	S4
## 13	B_D140314_SGSDSsample5_R01_MHRM		S5	S5
## 14	B_D140314_SGSDSsample5_R02_MHRM		S5	S5
## 15	B_D140314_SGSDSsample5_R03_MHRM		S5	S5
## 16	B_D140314_SGSDSsample6_R01_MHRM		S6	S6
## 17	B_D140314_SGSDSsample6_R02_MHRM		S6	S6
## 18	B_D140314_SGSDSsample6_R03_MHRM		S6	S6
## 19	B_D140314_SGSDSsample7_R01_MHRM		S7	S7
## 20	B_D140314_SGSDSsample7_R02_MHRM		S7	S7
## 21	B_D140314_SGSDSsample7_R03_MHRM		S7	S7
## 22	B_D140314_SGSDSsample8_R01_MHRM		S8	S8
## 23	B_D140314_SGSDSsample8_R02_MHRM		S8	S8
## 24	B_D140314_SGSDSsample8_R03_MHRM		S8	S8

5.2.2 Preprocessing with DIA experiment from Spectronaut output

The output from Spectronaut should look like below.

```
head(raw)
```

##	R.Condition	R.FileName	R.Replicate	PG.ProteinAccessions
## 1	SGSDSsample1 B_D140314_SGSDSsample1_R01_MHRM		1	A0A0B4J2A2
## 2	SGSDSsample1 B_D140314_SGSDSsample1_R01_MHRM		1	A0A0B4J2A2
## 3	SGSDSsample1 B_D140314_SGSDSsample1_R01_MHRM		1	A0A0B4J2A2
## 4	SGSDSsample1 B_D140314_SGSDSsample1_R01_MHRM		1	A0A0B4J2A2

```

## 5 SGSDSsample1 B_D140314_SGSDSsample1_R01_MHRM 1 AOA0B4J2A2
## 6 SGSDSsample1 B_D140314_SGSDSsample1_R01_MHRM 1 AOA0B4J2A2
## PG.ProteinGroups PG.Qvalue PG.Quantity PEP.GroupingKey PEP.StrippedSequence
## 1 AOA0B4J2A2 0 4662586 IIPGFMCQGGDFTR IIPGFMCQGGDFTR
## 2 AOA0B4J2A2 0 4662586 IIPGFMCQGGDFTR IIPGFMCQGGDFTR
## 3 AOA0B4J2A2 0 4662586 IIPGFMCQGGDFTR IIPGFMCQGGDFTR
## 4 AOA0B4J2A2 0 4662586 IIPGFMCQGGDFTR IIPGFMCQGGDFTR
## 5 AOA0B4J2A2 0 4662586 IIPGFMCQGGDFTR IIPGFMCQGGDFTR
## 6 AOA0B4J2A2 0 4662586 IIPGFMCQGGDFTR IIPGFMCQGGDFTR
## PEP.Quantity EG.iRTPredicted EG.Library
## 1 3303864 59.83870 S072_GSDS_DpD_Pulsar_Full.xls
## 2 3303864 59.83870 S072_GSDS_DpD_Pulsar_Full.xls
## 3 3303864 59.83870 S072_GSDS_DpD_Pulsar_Full.xls
## 4 3303864 59.83870 S072_GSDS_DpD_Pulsar_Full.xls
## 5 3303864 59.83870 S072_GSDS_DpD_Pulsar_Full.xls
## 6 3303864 45.20149 S072_GSDS_DpD_Pulsar_Full.xls
## EG.ModifiedSequence EG.PrecursorId
## 1 _IIPGFMC[+C2+H3+N+0]QGGDFTR_ _IIPGFMC[+C2+H3+N+0]QGGDFTR_.2
## 2 _IIPGFMC[+C2+H3+N+0]QGGDFTR_ _IIPGFMC[+C2+H3+N+0]QGGDFTR_.2
## 3 _IIPGFMC[+C2+H3+N+0]QGGDFTR_ _IIPGFMC[+C2+H3+N+0]QGGDFTR_.2
## 4 _IIPGFMC[+C2+H3+N+0]QGGDFTR_ _IIPGFMC[+C2+H3+N+0]QGGDFTR_.2
## 5 _IIPGFMC[+C2+H3+N+0]QGGDFTR_ _IIPGFMC[+C2+H3+N+0]QGGDFTR_.2
## 6 _IIPGFM[+0]C[+C2+H3+N+0]QGGDFTR_ _IIPGFM[+0]C[+C2+H3+N+0]QGGDFTR_.2
## EG.Qvalue FG.Charge
## 1 3.355320e-12 2
## 2 3.355320e-12 2
## 3 3.355320e-12 2
## 4 3.355320e-12 2
## 5 3.355320e-12 2
## 6 4.455933e-14 2
##
## 1 _IIPGFMC[+C2+H3+N+0]QGGDFTR_;AOA0B4J2A2.2;59.8387F:\\Data\\HRM_GS\\S072_GSDS_DpD_Pulsar_Full.
## 2 _IIPGFMC[+C2+H3+N+0]QGGDFTR_;AOA0B4J2A2.2;59.8387F:\\Data\\HRM_GS\\S072_GSDS_DpD_Pulsar_Full.
## 3 _IIPGFMC[+C2+H3+N+0]QGGDFTR_;AOA0B4J2A2.2;59.8387F:\\Data\\HRM_GS\\S072_GSDS_DpD_Pulsar_Full.
## 4 _IIPGFMC[+C2+H3+N+0]QGGDFTR_;AOA0B4J2A2.2;59.8387F:\\Data\\HRM_GS\\S072_GSDS_DpD_Pulsar_Full.
## 5 _IIPGFMC[+C2+H3+N+0]QGGDFTR_;AOA0B4J2A2.2;59.8387F:\\Data\\HRM_GS\\S072_GSDS_DpD_Pulsar_Full.
## 6 _IIPGFM[+0]C[+C2+H3+N+0]QGGDFTR_;AOA0B4J2A2.2;45.20149F:\\Data\\HRM_GS\\S072_GSDS_DpD_Pulsar_Full.
## FG.PrecMz FG.Quantity F.Charge F.FrgIon F.FrgLossType F.FrgMz F.FrgNum
## 1 799.8763 3027798.8 2 y12 no loss 686.7923 12
## 2 799.8763 3027798.8 1 y8 no loss 940.3942 8
## 3 799.8763 3027798.8 1 y6 no loss 652.3049 6
## 4 799.8763 3027798.8 1 y3 no loss 423.2350 3
## 5 799.8763 3027798.8 1 y7 no loss 780.3635 7
## 6 807.8738 257682.8 2 y12 no loss 694.7897 12
## F.FrgType F.ExcludedFromQuantification F.NormalizedPeakArea
## 1 y False 1880069.5
## 2 y False 382481.4
## 3 y False 331427.1
## 4 y False 244019.3
## 5 y False 189801.4
## 6 y False 197765.4
## F.NormalizedPeakHeight F.PeakArea F.PeakHeight
## 1 6854814.4 1687479.6 6152623.5
## 2 1320585.1 343300.9 1185307.5

```

## 3	1221639.0	297476.5	1096497.2
## 4	891761.8	219022.5	800411.9
## 5	673177.2	170358.6	604218.5
## 6	431701.5	174551.3	381027.4

The input data for MSstats is required to contain variables of ProteinName, PeptideSequence, PrecursorCharge, FragmentIon, ProductCharge, IsotopeLabelType, Condition, BioReplicate, Run, Intensity. These variable names should be fixed. Therefore, we need to get subset of useful columns and to rename them. Also several filtering steps are required. SpectronauttoMSstatsFormat function helps pre-processing for making right format of MSstats input from Spectronaut output. First, it uses only noLoss from F.FrgLossType. If not, multiple measurements for each feature and run can be happen. Spectronaut provides the column named F.ExcludedFromQuantification based on XIC quality such as interference between chromatographies. Only features with F.ExcludedFromQuantification == 'False' should be used. PG.ProteinGroups is used for ProteinName. EG.ModifiedSequence is used for PeptideSequence. FG.Charge is used for PrecursorCharge. F.FrgIon is used for FragmentIon. F.Charge is used for ProductCharge. F.PeakArea with default option is used for Intensity. Then several filtering steps will be performed.

Here is the summary of pre-processing steps for SWATH/DIA experiment in SpectronauttoMSstatsFormat function.

SpectronauttoMSstatsFormat

- Extract essential information(columns)
- Use only 'noLoss' in F.FrgLossType
- Filter by PG.Qvalue
- Remove shared peptides
- Aggregate multiple measurements per feature and run
- Remove features with all missing values or with less than 3 measurements across MS runs
- Remove protein with only one feature
- **Add annotation for experimental design** : Group, biological replicate, fraction information per MS run

```
## check options for converting format
?SpectronauttoMSstatsFormat
```

```
quant <- SpectronauttoMSstatsFormat(raw,
                                     annotation = annot,
                                     filter_with_Qvalue = TRUE, ## same as default
                                     qvalue_cutoff = 0.01, ## same as default
                                     removeProtein_with1Feature = TRUE)
```

```
## ** Intensities with great than 0.01 in PG.Qvalue are replaced with NA.
## ** Intensities with great than 0.01 in EG.Qvalue are replaced with zero.
## ** 175 features have all NAs or zero intensity values and are removed.
## ** All peptides are unique peptides in proteins.
## ** 38 features have 1 or 2 intensities across runs and are removed.
## ** All proteins have at least two features.
## ** No multiple measurements in a feature and a run.
```

This function shows the progress. The output of `SpectronauttoMSstatsFormat`, called `quant`, is ready for next step.

```
## now 'quant' is ready for MSstats
head(quant)
```

```
##      ProteinName      PeptideSequence PrecursorCharge FragmentIon
## 1  A0A0B4J2A2      _IIPGFMC[+C2+H3+N+O] QGGDFTR_      2          y12
## 2  A0A0B4J2A2      _IIPGFMC[+C2+H3+N+O] QGGDFTR_      2          y8
## 3  A0A0B4J2A2      _IIPGFMC[+C2+H3+N+O] QGGDFTR_      2          y6
## 4  A0A0B4J2A2      _IIPGFMC[+C2+H3+N+O] QGGDFTR_      2          y3
## 5  A0A0B4J2A2      _IIPGFMC[+C2+H3+N+O] QGGDFTR_      2          y7
## 6  A0A0B4J2A2 _IIPGFM[+O]C[+C2+H3+N+O] QGGDFTR_      2          y12
##      ProductCharge IsotopeLabelType Condition BioReplicate
## 1                2                  L      S1          S1
## 2                1                  L      S1          S1
## 3                1                  L      S1          S1
## 4                1                  L      S1          S1
## 5                1                  L      S1          S1
## 6                2                  L      S1          S1
##
##              Run Intensity
## 1 B_D140314_SGSDSsample1_R01_MHRM 1687479.6
## 2 B_D140314_SGSDSsample1_R01_MHRM 343300.9
## 3 B_D140314_SGSDSsample1_R01_MHRM 297476.5
## 4 B_D140314_SGSDSsample1_R01_MHRM 219022.5
## 5 B_D140314_SGSDSsample1_R01_MHRM 170358.6
## 6 B_D140314_SGSDSsample1_R01_MHRM 174551.3
```

5.2.3 Different options for Spectronaut output of DIA experiment in dataProcess

In `dataProcess`, users need to use `censoredInt='0'` for Spectronaut output. Spectronaut output generates very few number of NA. After applying Qvalue, zero intensities will be generated and those should be imputed.

```
spectronaut.proposed <- dataProcess(quant,
                                   normalization='equalizeMedian',
                                   summaryMethod="TMP",
                                   cutoffCensored="minFeature",
                                   censoredInt="0",
                                   MBimpute=TRUE,
                                   maxQuantileforCensored=0.999)
```

Further steps is the same as in general workflow (section 4.1).

5.3 Suggested workflow with DIA-Umpire output for DIA

This section describes steps and considerations to properly format data processed by DIA-Umpire for SWATH/DIA experiments, prior to the `MSstats` analysis. In the following example, the raw files for the benchmark dataset (Bruderer et al. 2015) are quantified by DIA-Umpire. The datasets and details for data processing are available in MassIVE.quant, MSV000081828, Reanalysis : RMSV000000252.3

5.3.1 Load DIA-Umpire output

We first load and access the dataset processed by DIA-Umpire.

```
# Read output from DIA_Umpire : three output from different levels.
raw.frag <- read.csv('dia_diaumpire/Bruderer2015_DIA_DIAumpire_input_FragSummary.xls', sep="\t")

raw.pep <- read.csv('dia_diaumpire/Bruderer2015_DIA_DIAumpire_input_PeptideSummary.xls', sep="\t")

raw.pro <- read.csv('dia_diaumpire/Bruderer2015_DIA_DIAumpire_input_ProtSummary.xls', sep="\t")
```

One file is for annotation information, required to fill in Condition and BioReplicate for corresponding Run information. Users have to prepare as csv or txt file like 'Bruderer2015_DIA_DIAumpire_annotation.csv', which includes Run, Condition, and BioReplicate information, and load it in R.

```
## Read in annotation including condition and biological replicates
annot <- read.csv("dia_diaumpire/Bruderer2015_DIA_DIAumpire_annotation.csv", header = TRUE)
annot
```

```
##      Condition BioReplicate      Run
## 1 SGSDSsample1 SGSDSsample1 B_D140314_SGSDSsample1_R01_MHRM_TO
## 2 SGSDSsample1 SGSDSsample1 B_D140314_SGSDSsample1_R02_MHRM_TO
## 3 SGSDSsample1 SGSDSsample1 B_D140314_SGSDSsample1_R03_MHRM_TO
## 4 SGSDSsample2 SGSDSsample2 B_D140314_SGSDSsample2_R01_MHRM_TO
## 5 SGSDSsample2 SGSDSsample2 B_D140314_SGSDSsample2_R02_MHRM_TO
## 6 SGSDSsample2 SGSDSsample2 B_D140314_SGSDSsample2_R03_MHRM_TO
## 7 SGSDSsample3 SGSDSsample3 B_D140314_SGSDSsample3_R01_MHRM_TO
## 8 SGSDSsample3 SGSDSsample3 B_D140314_SGSDSsample3_R02_MHRM_TO
## 9 SGSDSsample3 SGSDSsample3 B_D140314_SGSDSsample3_R03_MHRM_TO
## 10 SGSDSsample4 SGSDSsample4 B_D140314_SGSDSsample4_R01_MHRM_TO
## 11 SGSDSsample4 SGSDSsample4 B_D140314_SGSDSsample4_R02_MHRM_TO
## 12 SGSDSsample4 SGSDSsample4 B_D140314_SGSDSsample4_R03_MHRM_TO
## 13 SGSDSsample5 SGSDSsample5 B_D140314_SGSDSsample5_R01_MHRM_TO
## 14 SGSDSsample5 SGSDSsample5 B_D140314_SGSDSsample5_R02_MHRM_TO
## 15 SGSDSsample5 SGSDSsample5 B_D140314_SGSDSsample5_R03_MHRM_TO
## 16 SGSDSsample6 SGSDSsample6 B_D140314_SGSDSsample6_R01_MHRM_TO
## 17 SGSDSsample6 SGSDSsample6 B_D140314_SGSDSsample6_R02_MHRM_TO
## 18 SGSDSsample6 SGSDSsample6 B_D140314_SGSDSsample6_R03_MHRM_TO
## 19 SGSDSsample7 SGSDSsample7 B_D140314_SGSDSsample7_R01_MHRM_TO
## 20 SGSDSsample7 SGSDSsample7 B_D140314_SGSDSsample7_R02_MHRM_TO
## 21 SGSDSsample7 SGSDSsample7 B_D140314_SGSDSsample7_R03_MHRM_TO
## 22 SGSDSsample8 SGSDSsample8 B_D140314_SGSDSsample8_R01_MHRM_TO
## 23 SGSDSsample8 SGSDSsample8 B_D140314_SGSDSsample8_R02_MHRM_TO
## 24 SGSDSsample8 SGSDSsample8 B_D140314_SGSDSsample8_R03_MHRM_TO
```

5.3.2 Preprocessing with DIA experiment from DIA-Umpire output

Here is the summary of pre-processing steps for DIA experiment in DIAUmpiretoMSstatsFormat function.

DIAUmpiretoMSstatsFormat

- Get selected fragments from DIA-Umpire
- Get selected peptides from DIA-Umpire
- Subtract the peak intensities for selected peptides and fragments
- Change the data format
- Remove shared peptides
- Aggregate multiple measurements per feature and run
- Remove features with all missing values or with less than 3 measurements across MS runs
- Remove protein with only one feature
- **Add annotation for experimental design** : Group, biological replicate, fraction information per MS run

Options for DIAUmpiretoMSstatsFormat

- **raw.frag** : name of FragSummary_date.xls data, which includes feature-level data.
- **raw.pep** : name of PeptideSummary_date.xls data, which includes selected fragments information.
- **raw.pro** : name of ProteinSummary_date.xls data, which includes selected peptides information.
- **annotation** : name of 'annotation.txt' data which includes Raw.file, Condition, BioReplicate, and Run information.
- **useSelectedFrag** : TRUE (default) will use the selected fragment for each peptide. 'Selected_fragments' column is required.
- **useSelectedPep** : TRUE (default) will use the selected peptide for each protein. 'Selected_peptides' column is required.
- **summaryforMultipleRows** : max(default), sum, or mean. MSstats assumes that there is only one measurement (peak intensity) for one feature and one run. When there are multiple measurements for certain feature and certain run, MSstats need to know which measurements need to be used for further analysis. Users can use highest(max), sum or mean among multiple measurements for one feature and one run.
- **fewMeasurement** : remove or keep the feature with few measurements.
 - 'remove' : (default) remove the features that have 1 or 2 measurements across runs.
 - 'keep' : keep all the features. However, it could generate the error in the step for fitting the statistical model.
- **removeProtein_with1Feature** : TRUE will remove the proteins which have only 1 feature, which is the combination of peptide, precursor charge, fragment and charge. FALSE is default.

For further details, visit the help file using the following code.

```
## check options for converting format
?DIAUmpiretoMSstatsFormat
```

```
quant <- DIAUmpiretoMSstatsFormat(raw.frag, raw.pep, raw.pro,
                                  annot,
                                  useSelectedFrag = TRUE,
                                  useSelectedPep = FALSE,
                                  removeProtein_with1Feature = TRUE)
```

```
## ** Got the selected fragments.
```

```
## ** Extract the data from selected fragments and/or peptides.
## ** 0 features have all NAs or zero intensity values and are removed.
## ** 10824 features have 1 or 2 intensities across runs and are removed.
## ** 26 proteins, which have only one feature in a protein, are removed among 3645 proteins.
## ** No multiple measurements in a feature and a run.
```

This function shows the progress. The output of `DIAUmpiretoMSstatsFormat`, called `quant`, is ready for next step.

```
## now 'quant' is ready for MSstats
head(quant)
```

```
##                               Run                               ProteinName
## 1 B_D140314_SGSDSsample1_R01_MHRM_T0 sp|A0A0B4J2A2|PAL4C_HUMAN
## 2 B_D140314_SGSDSsample1_R01_MHRM_T0 sp|A0A0B4J2A2|PAL4C_HUMAN
## 3 B_D140314_SGSDSsample1_R01_MHRM_T0 sp|A0A0B4J2A2|PAL4C_HUMAN
## 4 B_D140314_SGSDSsample1_R01_MHRM_T0 sp|A0A0B4J2A2|PAL4C_HUMAN
## 5 B_D140314_SGSDSsample1_R01_MHRM_T0 sp|A0A0B4J2A2|PAL4C_HUMAN
## 6 B_D140314_SGSDSsample1_R01_MHRM_T0 sp|A0A0B4J2A2|PAL4C_HUMAN
##                               PeptideSequence FragmentIon Intensity
## 1 HTGSGILSMANAGPNTNGSQFFI[57.021(C)]CTAK_3      y10_1      NA
## 2 HTGSGILSMANAGPNTNGSQFFI[57.021(C)]CTAK_3      y11_1      NA
## 3 HTGSGILSMANAGPNTNGSQFFI[57.021(C)]CTAK_3       y4_1      NA
## 4 HTGSGILSMANAGPNTNGSQFFI[57.021(C)]CTAK_3       y5_1      NA
## 5 HTGSGILSMANAGPNTNGSQFFI[57.021(C)]CTAK_3       y6_1      NA
## 6 HTGSGILSMANAGPNTNGSQFFI[57.021(C)]CTAK_3       y7_1      NA
##   PrecursorCharge ProductCharge IsotopeLabelType   Condition BioReplicate
## 1              NA             NA                L SGSDSsample1 SGSDSsample1
## 2              NA             NA                L SGSDSsample1 SGSDSsample1
## 3              NA             NA                L SGSDSsample1 SGSDSsample1
## 4              NA             NA                L SGSDSsample1 SGSDSsample1
## 5              NA             NA                L SGSDSsample1 SGSDSsample1
## 6              NA             NA                L SGSDSsample1 SGSDSsample1
```

5.3.3 Different options for DIA-Umpire output of DIA experiment in `dataProcess`

In `dataProcess`, users need to use `censoredInt='NA'` for DIA-Umpire output.

```
diaumpire.proposed <- dataProcess(quant,
  normalization='equalizeMedian',
  summaryMethod="TMP",
  cutoffCensored="minFeature",
  censoredInt="NA", ## !! important
  MBimpute=TRUE,
  maxQuantileforCensored=0.999)
```

Further steps is the same as in general workflow (section 4.1).

5.4 Suggested workflow with OpenSWATH output for SWATH

This section describes steps and considerations to properly format data processed by OpenSWATH for SWATH experiments, prior to the MSstats analysis. In the following example, the dataset processed and quantified by OpenSWATH and available as supplementary in (Röst et al. 2014) is used. The datasets and details for data processing are available in MassIVE.quant, MSV000081829, Reanalysis : RMSV000000253.2

5.4.1 Load OpenSWATH output

```
## Read fragment-level data
```

```
raw <- read.csv("dia_openswath/Rost2014_DIA_OpenSWATH_input.txt", sep="\t")
head(raw)
```

```
##                                                                 transition_group
## 1      AQUA4SWATH_YeastB_GSMADVPK(UniMod:259)/2_run0_split_napedro_L120417_001_SW_combined.feature
## 2 AQUA4SWATH_HMLangeA_APIPTALDSTDSSK(UniMod:259)/2_run0_split_napedro_L120417_001_SW_combined.feature
## 3  AQUA4SWATH_HMLangeA_DITAFDETLFR(UniMod:267)/2_run0_split_napedro_L120417_001_SW_combined.feature
## 4  AQUA4SWATH_HMLangeA_LNTIYQNDLTK(UniMod:259)/2_run0_split_napedro_L120417_001_SW_combined.feature
## 5 AQUA4SWATH_HMLangeB_GDSSLLLVAVTEVK(UniMod:259)/2_run0_split_napedro_L120417_001_SW_combined.feature
## 6 AQUA4SWATH_HMLangeB_ITVDDSDQGANAK(UniMod:259)/2_run0_split_napedro_L120417_001_SW_combined.feature
##   decoy main_var_xx_swath_prelim_score var_bseries_score
## 1 FALSE                1.2878970                1
## 2 FALSE                0.1706103                5
## 3 FALSE                1.1339196                2
## 4 FALSE                1.7598289                2
## 5 FALSE               -0.4115174                2
## 6 FALSE                1.0441584                1
##   var_elution_model_fit_score var_intensity_score var_isotope_correlation_score
## 1                0.9342550                0.453363758                0.9621812
## 2                0.8913419                0.049083725                0.6132449
## 3                0.9999998                0.009851618                0.3554128
## 4                0.9159949                0.095405653                0.9065832
## 5                0.9574775                0.021453197                0.9215580
## 6                0.9746000                0.022167918                0.4878706
##   var_isotope_overlap_score var_library_corr var_library_rmsd var_log_sn_score
## 1                0.0000000        -0.1055161        0.27253090        2.6731795
## 2                0.06965649        -0.5652690        0.30030363        0.4186867
## 3                0.22701794        -0.6904735        0.07683261        0.7115797
## 4                0.10497624        0.9654186        0.21724993        2.2435553
## 5                0.02678354        -0.2470041        0.30010703        0.3742150
## 6                0.84251969        0.8373865        0.08390017        1.5353173
##   var_massdev_score var_massdev_score_weighted var_norm_rt_score
## 1         10.605493                8.722201        0.037829778
## 2          3.506636                1.525251        0.053104479
## 3          1.922673                1.770099        0.059778618
## 4          6.326571                7.671188        0.037207495
## 5         10.654442               13.899227        0.027236154
## 6          5.657834                3.023133        0.007901597
##   var_xcorr_coelution var_xcorr_coelution_weighted var_xcorr_shape
## 1          3.3763883                1.9254634        0.8278485
## 2          2.9055453                0.9019211        0.7211294
## 3          0.9163978                0.4743998        0.8000000
## 4          2.6757296                1.1959883        0.8199334
```

```

## 5          1.8944289          0.7742008          0.8040568
## 6          1.8595018          0.2477762          0.8121000
##   var_xcorr_shape_weighted var_yseries_score
## 1          0.7850388          2
## 2          0.7406344          3
## 3          0.7628001          3
## 4          0.7884498          2
## 5          0.8094922          0
## 6          0.8723227          0
##
##                                     transition_group
## 1          AQUA4SWATH_YeastB_GSMADVPK(UniMod:259)/2_run0_split_napedro_L120417_001_SW_combined.feature
## 2 AQUA4SWATH_HMLangeA_APIPTALDSTDSSK(UniMod:259)/2_run0_split_napedro_L120417_001_SW_combined.feature
## 3   AQUA4SWATH_HMLangeA_DITAFDETLFR(UniMod:267)/2_run0_split_napedro_L120417_001_SW_combined.feature
## 4   AQUA4SWATH_HMLangeA_LNTIYQNDLTK(UniMod:259)/2_run0_split_napedro_L120417_001_SW_combined.feature
## 5 AQUA4SWATH_HMLangeB_GDSSLLAVTEVK(UniMod:259)/2_run0_split_napedro_L120417_001_SW_combined.feature
## 6 AQUA4SWATH_HMLangeB_ITVDDSDQGANAK(UniMod:259)/2_run0_split_napedro_L120417_001_SW_combined.feature
##   run_id          filename          RT
## 1      0 split_napedro_L120417_001_SW_combined.featureXML 1024.232
## 2      0 split_napedro_L120417_001_SW_combined.featureXML 2525.123
## 3      0 split_napedro_L120417_001_SW_combined.featureXML 4923.724
## 4      0 split_napedro_L120417_001_SW_combined.featureXML 2396.661
## 5      0 split_napedro_L120417_001_SW_combined.featureXML 4217.872
## 6      0 split_napedro_L120417_001_SW_combined.featureXML 1058.414
##           id          Sequence          FullPeptideName Charge          m.z
## 1 f_1353009549277083696      GSMADVPK      GSMADVPK(UniMod:259)          2 406.707
## 2 f_17169785622655779335 APIPTALDSTDSSK APIPTALDSTDSSK(UniMod:259)          2 662.348
## 3 f_14843615568932246264   DITAFDETLFR   DITAFDETLFR(UniMod:267)          2 669.334
## 4 f_2705275134670444755   LNTIYQNDLTK   LNTIYQNDLTK(UniMod:259)          2 665.858
## 5 f_9381457823485960609   GDSSLLAVTEVK   GDSSLLAVTEVK(UniMod:259)          2 670.382
## 6 f_4286219334306522917   ITVDDSDQGANAK   ITVDDSDQGANAK(UniMod:259)          2 671.322
##   Intensity          ProteinName assay_rt  delta_rt leftWidth  norm_RT
## 1    228484   AQUA4SWATH_YeastB 1160.081 -135.84854   1003.92 -19.18298
## 2    11528   AQUA4SWATH_HMLangeA 2343.868  181.25584   2513.89  24.41045
## 3    1784   AQUA4SWATH_HMLangeA 4711.441  212.28267   4920.69  94.07786
## 4    41457   AQUA4SWATH_HMLangeA 2525.725 -129.06375   2384.16  20.67925
## 5    4107   AQUA4SWATH_HMLangeB 4306.552  -88.67945   4214.01  73.57638
## 6     762   AQUA4SWATH_HMLangeB 1091.456  -33.04119   1045.90 -18.19016
##   nr_peaks peak_apices_sum rightWidth  rt_score  sn_ratio total_xic
## 1      4          23302    1061.95 3.7829778 14.485954   503975
## 2      4          1673    2537.79 5.3104479  1.519964   234864
## 3      4           902    4924.10 5.9778618  2.037207   181087
## 4      4          8532    2411.47 3.7207495  9.426787   434534
## 5      4          1065    4227.67 2.7236154  1.453850   191440
## 6      4           291    1062.97 0.7901597  4.642799    34374
##   dotprod_score library_dotprod library_manhattan manhatt_score
## 1    0.7223189    0.7550611    0.6275367    0.7230707
## 2    0.6545768    0.7396598    0.6218496    0.8495160
## 3    0.7349210    0.9816014    0.1471980    0.7705343
## 4    0.7527736    0.8719077    0.5787663    0.7852044
## 5    0.6493268    0.7414075    0.6934909    0.8909054
## 6    0.4477453    0.9752496    0.2175255    1.0719205
##   xx_lda_prelim_score xx_swath_prelim_score
## 1          3.488009          0
## 2          1.007469          0

```


## 21	split_napedro_L120420_001_SW_combined.featureXML	512	512
## 22	split_napedro_L120420_002_SW_combined.featureXML	256	256
## 23	split_napedro_L120420_003_SW_combined.featureXML	128	128
## 24	split_napedro_L120420_004_SW_combined.featureXML	64	64
## 25	split_napedro_L120420_005_SW_combined.featureXML	32	32
## 26	split_napedro_L120420_006_SW_combined.featureXML	16	16
## 27	split_napedro_L120420_007_SW_combined.featureXML	8	8
## 28	split_napedro_L120420_008_SW_combined.featureXML	4	4
## 29	split_napedro_L120420_009_SW_combined.featureXML	2	2
## 30	split_napedro_L120420_010_SW_combined.featureXML	1	1
##	Run		
## 1	split_napedro_L120417_001_SW_combined.featureXML		
## 2	split_napedro_L120417_002_SW_combined.featureXML		
## 3	split_napedro_L120417_003_SW_combined.featureXML		
## 4	split_napedro_L120417_004_SW_combined.featureXML		
## 5	split_napedro_L120417_005_SW_combined.featureXML		
## 6	split_napedro_L120417_006_SW_combined.featureXML		
## 7	split_napedro_L120417_007_SW_combined.featureXML		
## 8	split_napedro_L120417_008_SW_combined.featureXML		
## 9	split_napedro_L120417_009_SW_combined.featureXML		
## 10	split_napedro_L120417_010_SW_combined.featureXML		
## 11	split_napedro_L120419_001_SW_combined.featureXML		
## 12	split_napedro_L120419_002_SW_combined.featureXML		
## 13	split_napedro_L120419_003_SW_combined.featureXML		
## 14	split_napedro_L120419_004_SW_combined.featureXML		
## 15	split_napedro_L120419_005_SW_combined.featureXML		
## 16	split_napedro_L120419_006_SW_combined.featureXML		
## 17	split_napedro_L120419_007_SW_combined.featureXML		
## 18	split_napedro_L120419_008_SW_combined.featureXML		
## 19	split_napedro_L120419_009_SW_combined.featureXML		
## 20	split_napedro_L120419_010_SW_combined.featureXML		
## 21	split_napedro_L120420_001_SW_combined.featureXML		
## 22	split_napedro_L120420_002_SW_combined.featureXML		
## 23	split_napedro_L120420_003_SW_combined.featureXML		
## 24	split_napedro_L120420_004_SW_combined.featureXML		
## 25	split_napedro_L120420_005_SW_combined.featureXML		
## 26	split_napedro_L120420_006_SW_combined.featureXML		
## 27	split_napedro_L120420_007_SW_combined.featureXML		
## 28	split_napedro_L120420_008_SW_combined.featureXML		
## 29	split_napedro_L120420_009_SW_combined.featureXML		
## 30	split_napedro_L120420_010_SW_combined.featureXML		

5.4.2 Preprocessing with DIA experiment from OpenSWATH output

The output from OpenSWATH should look like below.

```
head(raw)
```

Here is the summary of pre-processing steps for SWATH/DIA experiment in `OpenSWATHtoMSstatsFormat` function.

OpenSWATHtoMSstatsFormat

- Remove the decoys
- Filter by `m_score`
- Change data format
- Remove shared peptides
- Aggregate multiple measurements per feature and run : max or mean
- Remove features with all missing values or with less than 3 measurements across MS runs
- Remove protein with only one feature
- **Add annotation for experimental design** : Group, biological replicate, fraction information per MS run

Options for OpenSWATHtoMSstatsFormat

- **input** : name of MSstats input report from OpenSWATH, which includes feature-level data.
- **annotation** : name of 'annotation.txt' data which includes Raw.file, Condition, BioReplicate, and Run information. Run should be the same as 'filename'.
- **filter_with_mscore** : TRUE (default) will filter out the features that have greater than `mscore_cutoff` in `m_score` column. Those features will be removed.
- **mscore_cutoff** : Cutoff for `m_score`. default is 0.01.
- **useUniquePeptide** : TRUE (default) removes peptides that are assigned for more than one proteins. We assume to use unique peptide for each protein.
- **summaryforMultipleRows** : max(default), sum, or mean. MSstats assumes that there is only one measurement (peak intensity) for one feature and one run. When there are multiple measurements for certain feature and certain run, MSstats need to know which measurements need to be used for further analysis. Users can use highest(max), sum or mean among multiple measurements for one feature and one run.
- **fewMeasurement** : remove or keep the feature with few measurements.
 - 'remove' : (default) remove the features that have 1 or 2 measurements across runs.
 - 'keep' : keep all the features. However, it could generate the error in the step for fitting the statistical model.
- **removeProtein_with1Feature** : TRUE will remove the proteins which have only 1 peptide and charge. FALSE is default.

For further details, visit the help file using the following code.

```
## check options for converting format
?OpenSWATHtoMSstatsFormat
```

Now, we use `OpenSWATHtoMSstatsFormat` function for this example dataset. We chose to remove the proteins with only 1 feature.

```
quant <- OpenSWATHtoMSstatsFormat(raw,
                                   annotation = annot,
                                   filter_with_mscore = TRUE, ## same as default
                                   mscore_cutoff = 0.01, ## same as default
                                   removeProtein_with1Feature = TRUE)
```

```
## ** Features with great than 0.01 in m_score are removed.
```



```
## ** 0 features have all NAs or zero intensity values and are removed.
## ** All peptides are unique peptides in proteins.
## ** 31 features have 1 or 2 intensities across runs and are removed.
## ** All proteins have at least two features.
## ** No multiple measurements in a feature and a run.
```

This function shows the progress. The output of `OpenMStoMSstatsFormat`, called `quant`, is ready for next step.

```
## now 'quant' is ready for MSstats
head(quant)
```

```
##           ProteinName           PeptideSequence PrecursorCharge
## 1 AQUA4SWATH_HMLangeA ADSTGTLVITDPTR(UniMod_267)             2
## 2 AQUA4SWATH_HMLangeA ADSTGTLVITDPTR(UniMod_267)             2
## 3 AQUA4SWATH_HMLangeA ADSTGTLVITDPTR(UniMod_267)             2
## 4 AQUA4SWATH_HMLangeA ADSTGTLVITDPTR(UniMod_267)             2
## 5 AQUA4SWATH_HMLangeA ALGYEDATQALGR(UniMod_267)             2
## 6 AQUA4SWATH_HMLangeA ALGYEDATQALGR(UniMod_267)             2
##                                     FragmentIon ProductCharge
## 1 AQUA4SWATH_HMLangeA_ADSTGTLVITDPTR(UniMod_267)/2_y10      NA
## 2 AQUA4SWATH_HMLangeA_ADSTGTLVITDPTR(UniMod_267)/2_y7      NA
## 3 AQUA4SWATH_HMLangeA_ADSTGTLVITDPTR(UniMod_267)/2_y8      NA
## 4 AQUA4SWATH_HMLangeA_ADSTGTLVITDPTR(UniMod_267)/2_y9      NA
## 5 AQUA4SWATH_HMLangeA_ALGYEDATQALGR(UniMod_267)/2_y7      NA
## 6 AQUA4SWATH_HMLangeA_ALGYEDATQALGR(UniMod_267)/2_y8      NA
## IsotopeLabelType Condition BioReplicate
## 1                L        512         512
## 2                L        512         512
## 3                L        512         512
## 4                L        512         512
## 5                L        512         512
## 6                L        512         512
##                                     Run Intensity
## 1 split_napedro_L120417_001_SW_combined.featureXML          0
## 2 split_napedro_L120417_001_SW_combined.featureXML          0
## 3 split_napedro_L120417_001_SW_combined.featureXML          0
## 4 split_napedro_L120417_001_SW_combined.featureXML          0
## 5 split_napedro_L120417_001_SW_combined.featureXML          0
## 6 split_napedro_L120417_001_SW_combined.featureXML          0
```

5.4.3 Different options for OpenSWATH output of DIA experiment in `dataProcess`

In `dataProcess`, users need to use `censoredInt='0'` for OpenSWATH output.

```
goldstandard.proposed <- dataProcess(quant,
                                     normalization='equalizeMedian',
                                     summaryMethod="TMP",
                                     cutoffCensored="minFeature",
                                     censoredInt="0",
                                     MBimpute=TRUE,
                                     maxQuantileforCensored=0.999)
```


Further steps is the same as in general workflow (section 4.1).

6. SRM analysis with MSstats

6.1 Suggested workflow for SRM

This section describes a typical workflow for SRM experiments with heavy labeled-isotope peptides. The example dataset, `SRMRawData` in `MSstats` is used for demonstration.

6.1.1 Preparing the data for MSstats input

The first step in using the `MSstats` is to format the data as described in Section 2. `SRMRawData` is already formatted for `MSstats` input.

```
# Check the first 6 rows in SRMRawData
head(SRMRawData)
```

```
##      ProteinName PeptideSequence PrecursorCharge FragmentIon ProductCharge
## 243      IDHC     ATDVIVPEEGELR             2          y7           NA
## 244      IDHC     ATDVIVPEEGELR             2          y7           NA
## 245      IDHC     ATDVIVPEEGELR             2          y8           NA
## 246      IDHC     ATDVIVPEEGELR             2          y8           NA
## 247      IDHC     ATDVIVPEEGELR             2          y9           NA
## 248      IDHC     ATDVIVPEEGELR             2          y9           NA
##      IsotopeLabelType Condition BioReplicate Run      Intensity
## 243                H         1      ReplA    1 84361.08350
## 244                L         1      ReplA    1  215.13526
## 245                H         1      ReplA    1 29778.10188
## 246                L         1      ReplA    1   98.02134
## 247                H         1      ReplA    1 17921.29255
## 248                L         1      ReplA    1   60.47029
```

6.1.2 Processing the data

It is the same workflow as described in section 4.1.2. Only difference is the normalization with heavy labeled isotope peptides.


Different normalization option Let's see the different normalization effect with SRM dataset including two proteins.

```
unique(SRMRawData$ProteinName)
```

```
## [1] IDHC PMG2
## 45 Levels: ACEA ACH1 ACON ADH1 ADH2 ADH4 ALDH6 ALF CISO1 CISO2 DHSA ... SUCB
```


No normalization No normalization is performed. If you had your own normalization before `MSstats`, you should use like below.

```
srm.nonorm <- dataProcess(SRMRawData, normalization=FALSE)
dataProcessPlots(srm.nonorm, type='QCplot', address='srm_noNorm_')
```

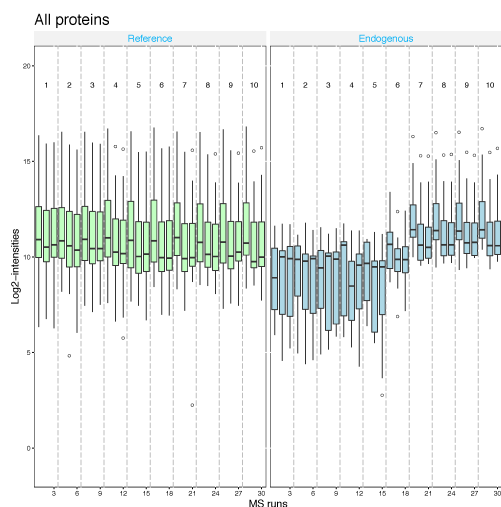


Equalize medians normalization The default option for normalization is ‘`equalizeMedians`’, where all the intensities in a run are shifted by a constant, to equalize the median of intensities across runs for label-free experiment. This normalization method is appropriate **when we can assume that the majority of proteins do not change across runs**. Be cautious when using the `equalizeMedians` option for a label-free dataset with only a small number of proteins. For label based experiment, `equalizeMedians` equalizes the median of reference intensities across runs and is generally proper even for a dataset with a small number of proteins.

```
srm.equalmed <- dataProcess(SRMRawData, normalization = 'equalizeMedians')
dataProcessPlots(srm.equalmed, type='QCplot', address='srm_equalM_')
```

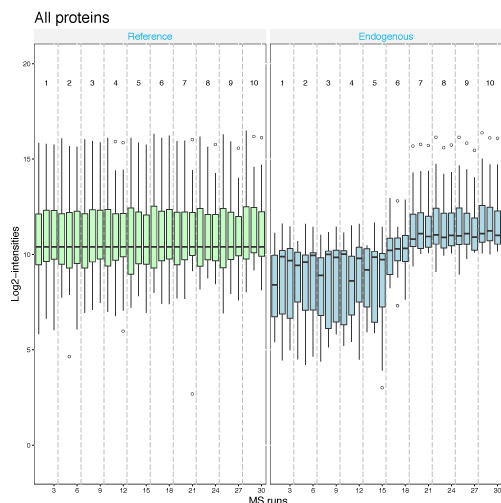

Quantile normalization The distribution of all the intensities in each run will become the same across runs for label-free experiment. For label-based experiment, the distribution of all the reference intensities will be become the same across runs and all the endogenous intensities are shifted by a constant corresponding to reference intensities.


```
srm.quantile <- dataProcess(SRMRawData, normalization='quantile')
dataProcessPlots(srm.quantile, type='QCplot', address='srm_quantile_')
```

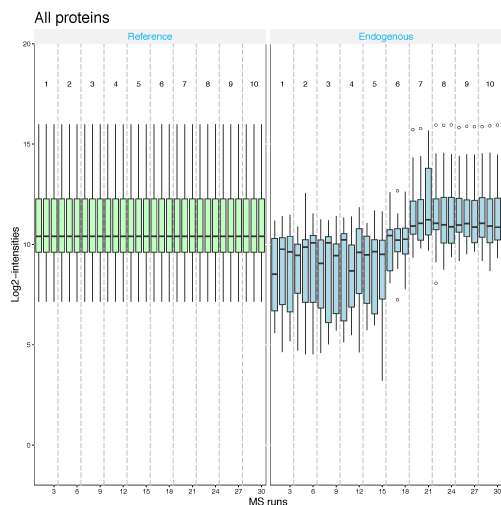

Global standards normalization : example 1 If you have a spiked in standard across all MS runs, you may set this to `globalStandards` and define the standard with `nameStandards` option. Global standard peptide or Protein names, which you can assume that they have the same abundance across MS runs, should be assigned in the vector for this option.


First, let's assume that PMG2 proteins is the spike-in protein and should be equal amount across MS runs.

```
srm.global.pmg2 <- dataProcess(SRMRawData, normalization = 'globalStandards',
                               nameStandards = 'PMG2')
dataProcessPlots(srm.global.pmg2, type='QCplot', address='srm_global_PMG2_')
```

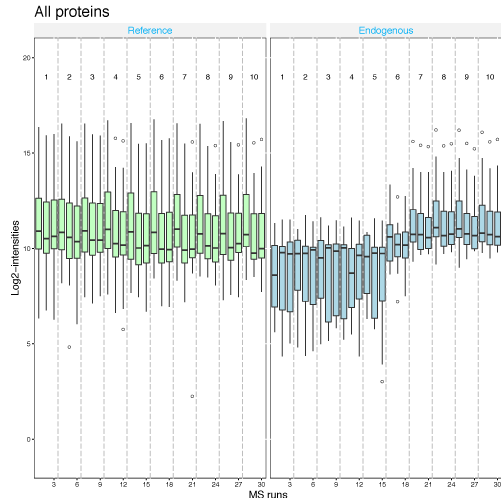

Second, let's assume that IDHC proteins is the spike-in protein and should be equal amount across MS runs.

```
srm.global.idhc <- dataProcess(SRMRawData, normalization = 'globalStandards',
                               nameStandards = 'IDHC')
dataProcessPlots(srm.global.idhc, type='QCplot', address='srm_global_IDHC_')
```

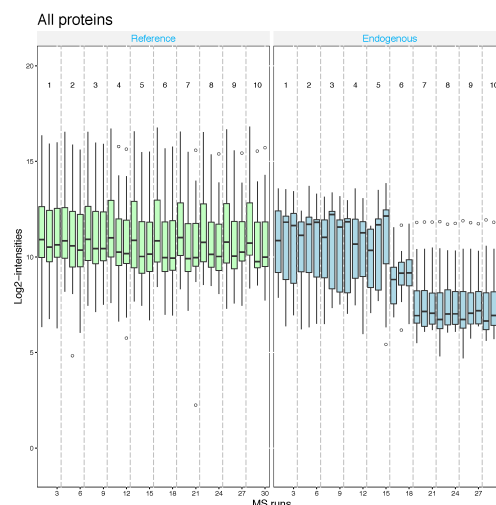

Global standards normalization : example 2

Further steps is the same as in general workflow (section 4.1).

Reference

Benjamini, Y., and Y. Hochberg. 1955. "Controlling the false discovery rate: a practical and powerful approach to multiple testing." *J.R. Statist. Soc. B* 57 (1): 289–300.

Bruderer, R., O. M. Bernhardt, T. Gandhi, S. M. Miladinović, L.-Y. Cheng, S. Messner, T. Ehrenberger, et al. 2015. "Extending the limits of quantitative proteome profiling with Data-Independent Acquisition and application to acetaminophen-treated three-dimensional liver microtissues." *Mol. Cell. Proteomics* 14 (5): 1400–1410.

Choi, M. and Eren-Dogu, Z. F. and Colangelo, C. and Cottrell, J. and Hoopmann, M. R. and Kapp, E. A. and Kim, S. and Lam, H. and Neubert, T. A. and Palmblad, M. and Phinney, B. S. and Weintraub, S. T. and MacLean, B. and Vitek, O. 2017. "ABRF Proteome Informatics Research Group (iPRG) 2015 Study: Detection of differentially abundant proteins in label-free quantitative LC-MS/MS experiments." *J. Proteome Res.* <https://doi.org/10.1021/acs.jproteome.6b00881>.

Cox, Jürgen, and Matthias Mann. 2008. "MaxQuant enables high peptide identification rates, individualized p.p.b.-range mass accuracies and proteome-wide protein quantification." *Nature Biotechnology* 26 (12): 1367–72.

MacLean, B., D. M. Tomazela, N. Shulman, M. Chambers, G. Finney, B. Frewen, R. Kern, D. L. Tabb, D. C. Liebner, and M. J. MacCoss. 2010. "Skyline: An open source document editor for creating and analyzing targeted proteomics experiments." *Bioinformatics* 26–27: 966.

Mueller, L. N., O. Rinner, A. Schmidt, S. Letarte, B. Bodenmiller, M.-Y. Brusniak, O. Vitek, R. Aebersold, and M. Müller. 2007. "SuperHirn - a novel tool for high resolution LC-MS-based peptide/protein profiling." *Proteomics* 7: 3470–80.

Röst, H. L., G. Rosenberger, P. Navarro, L. Gillet, S. M. Miladinović, O. T. Schubert, W. Wolski, et al. 2014. "OpenSWATH enables automated, targeted analysis of data-independent acquisition MS data." *Nat. Biotechnol.* 32 (3): 219–23.

Sturm, Marc, Andreas Bertsch, Clemens Gröpl, Andreas Hildebrandt, Rene Hussong, Eva Lange, Nico Pfeifer, et al. 2008. "OpenMS – An open-source software framework for mass spectrometry." *BMC Bioinformatics* 9 (1): 163.

Tsou, C.-C., D. Avtonomov, B. Larsen, M. Tucholska, H. Choi, A.-C. Gingras, and A. I. Nesvizhskii. 2015. “DIA-Umpire: comprehensive computational framework for data-independent acquisition proteomics.” *Nat. Methods* 12 (3): 258–64.